

Tetracycline assay method.

## FIELD OF THE INVENTION

This invention relates to a method for the determination of a tetracycline in a sample. The invention also concerns recombinant prokaryotic cells capable of emitting light in response to the existence of a tetracycline in a sample. Furthermore, the invention relates to novel DNA vectors useful for the construction of said prokaryotic cells.

## 10 BACKGROUND OF THE INVENTION

The publications and other materials used herein to illuminate the background of the invention, and in particular, cases to provide additional details respecting the practice, are incorporated by reference.

15 Whole cells can be used in methods based on the use of living cells or organisms as sensor tools of detection. Many of these methods utilize bacterial or yeast cells. Prokaryotic organisms and especially *Escherichia coli* bacterium are very well characterized and maps of genes and their sequences at nucleotide level are known. Therefore the behavior of the whole cell sensor can be better understood. Because  
20 of this fact it is also possible to develop analyte or group specific sensors utilizing different regulatory regions of genomes and also various microbial strains.

Whole cells can be utilized in biosensors which are devices consisting of 1) a sensor, 2) a recording unit and 3) a possible connector such as fiber optic guide  
25 between 1 and 2. The recording unit has several choices of what is the physical background of the measurement. It can be change in heat, conductance, color reaction, changes in fluorescent properties, emission of endogenous light from the sensor cells etc.

004240 29662560

- Antibiotics used as medicines against microbial invasion are detected from body fluids in order to study the dosage and penetration of the medicine. Often the effective therapeutic range of the antibiotic is rather narrow and the risks of overdosage might be too big. It is also important to measure the presence or concentration of antibiotics from meat and milk due to syndrome of allergic people. In the course of cheese production milk used as starting material should not contain antibiotics due to the fact that cheesemaking bacteria are not able to work on contaminated milk.
- Conventional tests for the measurement of toxic substances such as antimicrobial agents (antibiotics) are based on the inhibition of growth. Growth inhibition can be followed by monitoring the zone where the growth of microbes is inhibited on a nutrient agar plate around a disk onto which an antibiotic dilution was pipetted.
- Typical examples of agar diffusion tests are cylindrical, hole or disk methods. The difference in these tests is only restricted in the way the sample is applied on the agar and also the way the bacteria in the test is used. Another means is to follow the metabolism of the test organisms by estimating the intensity of a color reaction which is affected by the inhibitory antibiotic present and comparing it to the uninhibited control (e.g. the commercial products: Delvo Test<sup>TM</sup>, Brilliant black-reduction test, Charm Farm Test, Charm AIM-96 and Valio T101-test). Since microbiological methods utilize bacteria or their spores it is the sensitivity of the test bacteria which is of utmost importance. Thus far one had to make compromises in the choice of a suitable test strain since great sensitivity against antimicrobial agents and other characteristics needed for the test strain have not been common features for the same strain of bacteria. A major drawback when using microbes in antibiotic residue tests is slow and unsensitive performance. Since in these methods one always controls in a way or other the growth of the tester strain one cannot imagine

the test to be performed in an hour. This is due to the fact that the growth of the microbe is a slow phenomenon even at its fastest mode. Also in many cases microbes are in spores or freeze-dried, the regeneration of which makes the tests even more slow to perform.

5

## OBJECT AND SUMMARY OF THE INVENTION

The object of the invention is to provide a novel method of determining a tetracycline in a sample where said method is rapid and selective for tetracyclines, i.e. the method is able to distinguish tetracyclines from other antimicrobial agents.

10

According to one aspect of the invention a method for the determination of a tetracycline in a sample is provided, wherein the method is characterized in that

- the sample is brought into contact with prokaryotic cells encompassing a DNA vector including a nucleotide sequence encoding a light producing enzyme under transcriptional control of a tetracycline repressor and a tetracycline promoter,
  - detecting the luminescence emitted from the cells, and
  - comparing the emitted luminescence to the luminescence emitted from cells in a control containing no tetracycline
  - wherein a detectable luminescence higher than a luminescence of the control
- 20 indicates the presence of tetracycline in the sample.

According to another aspect, the invention concerns a recombinant prokaryotic cell which encompasses a DNA vector including a nucleotide sequence encoding a light producing enzyme, tetracycline repressor and tetracycline promoter.

25

According to yet another aspect, the invention concerns a plasmid which comprises either

004240 042400 004240 042400

- the luxCDABE genes (SEQ ID NO: 3), tetracycline repressor (TetR) (SEQ ID NO: 11) and tetracycline promotor (TetA) (SEQ ID NO: 9) from *Tn10*, or
- the insect luciferase gene (SEQ ID NO: 1), tetracycline repressor (TetR) (SEQ ID NO: 11) and tetracycline promotor (TetA) (SEQ ID NO: 9) from *Tn10*.

5

### BRIEF DESCRIPTION OF THE DRAWINGS

Figure 1a shows schematically the method according to this invention, where cells cloned with the plasmid pTetLux1 (SEQ ID NO: 3) are used.

- 10 Figure 1b shows schematically the method according to this invention, where cells cloned with the plasmid pTetLuc1 (SEQ ID NO: 1) are used.

Figure 1c shows schematically the production of the luciferase enzyme,

- 15 Figure 2 shows the plasmid pTetLux1 (SEQ ID NO: 3).

Figure 3 shows the plasmid pTetLuc1 (SEQ ID NO: 1).

- 20 Figure 4a shows the production of light (induction factor) versus concentration of tetracycline in samples for three different tetracyclines,

Figure 4b shows the production of light (induction factor) versus concentration of tetracycline in samples for further four different tetracyclines.

- 25 Figure 5 shows the effect of magnesium ions on the sensitivity of the method according to the invention.

Figure 6 illustrates possibilities of changing the assay window for the method of the invention by adjusting magnesium ion concentration and pH.

Figure 7 shows the induction factor versus tetracycline concentration when using freeze-dried *E. coli* in the determination of tetracycline.

Figure 8 shows a comparison of the assays based on using cells with the plasmid pTetLuc1 (SEQ ID NO: 1) and with the plasmid pTetLux1 (SEQ ID NO: 3).


Figure 9 shows induction factors versus antibiotic concentrations of a pig serum sample (cells *E. coli* K12, pTetLux1).

Figure 10 shows the effect of EDTA in a milk sample assay, and

Figure 11 shows the light emission versus time for an assay according to the invention.

## DETAILED DESCRIPTION OF THE INVENTION

The term "tetracycline" shall be understood to include any compound covered by the general structure formula

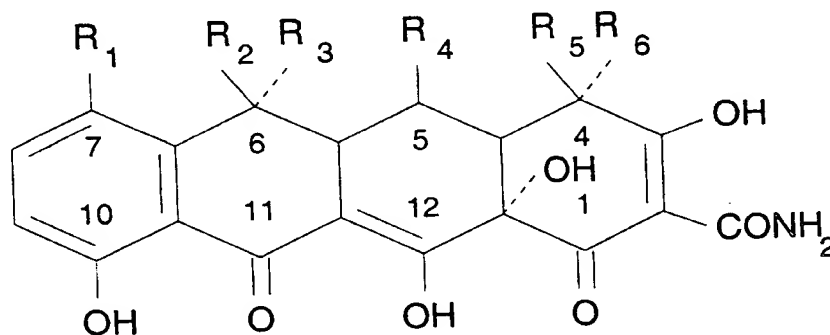

and particularly the specific commercially available compounds listed in the table below.

GENERIC NAME	R <sub>1</sub>	R <sub>2</sub>	R <sub>3</sub>	R <sub>4</sub>	R <sub>5</sub>	R <sub>6</sub>
Chlorotetracycline	Cl	OH	CH <sub>3</sub>	H	H	N(CH <sub>3</sub> ) <sub>2</sub>
Demethylchlorotetracycline	Cl	OH	H	H	H	N(CH <sub>3</sub> ) <sub>2</sub>
Doxycycline	H	H	CH <sub>3</sub>	OH	H	N(CH <sub>3</sub> ) <sub>2</sub>
Methacycline	H	CH <sub>3</sub>	H	OH	H	N(CH <sub>3</sub> ) <sub>2</sub>
Minocycline	N(CH <sub>3</sub> ) <sub>2</sub>	H	H	H	H	N(CH <sub>3</sub> ) <sub>2</sub>
Oxytetracycline	H	OH	CH <sub>3</sub>	OH	H	N(CH <sub>3</sub> ) <sub>2</sub>
Tetracycline	H	OH	CH <sub>3</sub>	H	H	N(CH <sub>3</sub> ) <sub>2</sub>

Furthermore, the term "tetracycline" shall be understood to cover the metabolic and other reformulation/decomposition products thereof.

5

The cells useful in the method of the invention are preferably *Escherichia coli*, which are stored in dried form, e.g. in lyophilized form before their use in the method according to the invention. Also freshly cultivated cells can be used.

- 10 According to a preferred embodiment, the DNA vector including a nucleotide sequence encoding a light producing enzyme is a plasmid containing the luxCDABE genes (SEQ ID NO: 3), tetracycline repressor (TetR) (SEQ ID NO: 11) and tetracycline promotor (TetA) (SEQ ID NO: 9) from transposon *Tn10*. Particularly preferable is the plasmid pTetLux1 (SEQ ID NO: 3).

15

According to another preferred embodiment, the DNA vector including a nucleotide sequence encoding a light producing enzyme is a plasmid containing the insect

luciferase gene, tetracycline repressor (TetR) (SEQ ID NO: 11) and tetracycline promotor (TetA) (SEQ ID NO: 9) from *Tn10*. In this case the substrate for insect luciferase reaction, D-luciferin, is added to the mixture of the sample and the cells in order to initiate the luminescence of the cells. The plasmid is preferably pTetLuc1

5 (SEQ ID NO: 1).

The method according to this invention is useful for the determination of tetracycline in various kinds of samples. As examples can be mentioned milk, fish, meat, infant formula, eggs, honey, vegetables, serum, plasma, whole blood or the

10 like.

sub  
ci } The luminescence of the cells is preferably measured using an X-ray or polaroid film, a CCD-camera (Charge Coupled Device), a liquid scintillation counter or, most preferably, a luminometer.

15

The sensitivity of this analysis method with respect to the tetracycline can be controlled by increasing or decreasing the concentration of divalent metal ions, e.g. magnesium ions, in the mixture of the sample and the cells, by adjusting the pH or by combined adjusting of the divalent metal ion concentration and the pH.

20 Increasing concentration of magnesium ions decreases the sensitivity and vice versa. Increasing pH will also cause a decreasing sensitivity. The sensitivity of the analysis with respect to the tetracycline can be increased by the use of cells which are especially antibiotic sensitive mutant strains. Chelating agents such as EDTA can be added to further sensitize the sensor system for tetracyclines.

25

Figures 1 show a schematic representation of a method based on specific detection of the presence of tetracyclines using microbial cells cloned with either the plasmid pTetLux1 (SEQ ID NO: 3) (Figure 1a) or with the plasmid pTetLuc1 (SEQ ID

004240 29662560

NO: 1) (Figure 1b). The figures show that cells containing either of the plasmids can be triggered to produce light by adding a chemical agent (a tetracycline). Light production is a consequence of tetracycline responsive promoter activation due to removal of the tet-repressor protein (SEQ ID NO: 11) leading to the production of  
5 luciferase specific mRNA and luciferase protein (SEQ ID NO: 2, 4-8) itself. The principle is demonstrated in Figure 1c. In case of the usage of full length bacterial luciferase operon (SEQ ID NO: 3) containing *luxC*, *luxD*, *luxA*, *luxB* and *luxE* genes (SEQ ID NO: 3) (Figure 1a), one is able to get light emission without addition of any substance. In case of insect (e.g. firefly) luciferase (SEQ ID NO: 2) (Figure  
10 1b), light is emitted only after the addition of D-luciferin. It should be noticed that the triggering of luciferase synthesis and light production commences immediately when the cells are introduced to the inducer molecules (tetracyclines). Therefore there is no need to use dividing cells and hence there is no need to use long cultivation of microbial cells such as the case is with conventional methods.  
15 Therefore, if needed, one can get results in minutes rather than in hours or days which is the case when conventional methods are used.

Figure 2a shows the plasmid pTetLux1 (SEQ ID NO: 3), in which the production of bacterial luciferase (SEQ ID NO: 4-8) of *Photorhabdus luminescens* (formerly  
20 *Xenorhabdus luminescens*; the lux-operon structure and the full-length nucleotide sequence of *P. luminescens* was published in Szittner, R. and Meighen, E. (1990) J. Biol. Chem. 265, 16581-16587) can be switched on by the addition of a chemical agent belonging to the tetracycline family of antimicrobial agents in a cloned *E. coli* bacterium. SEQ ID NO: 3 shows the nucleotide sequence of the plasmid pTetLux1.  
25 This plasmid construct is devised to contain the five genes from *P. luminescens* luciferase operon necessary for the light production without any additions of substrates, i.e. cells cloned with such a construct produce substrates endogenously. By incubating *E. coli* cells containing this plasmid (or any other microbial strain



where to similar regulation/reporter gene system is incorporated containing the necessary secondary regulatory sequences in the constructs such as correct ribosome binding region, transcriptional termination etc.) in the presence of very small amounts of tetracyclines one is able to obtain light production the intensity of which is proportional to the concentration of tetracycline used.

Any *E. coli* mutant strain and especially those strains having a mutation in the export/import machinery of the membranes or otherwise leaky character making it possible for large molecules to easily penetrate inside the cell would be beneficial to use in the method described in this invention. Also other gram-negative bacteria such as strains belonging to genus *Salmonella*, *Shigella*, *Enterobacter*, *Citrobacter*, *Klebsiella*, *Erwinia*, *Pseudomonas*, *Serratia* as well as gram-positive organisms such as those belonging to genus *Bacillus* (especially *B. subtilis*, *B. licheniformis*, *B. pumilus*, *B. globigii*, *B. natto*, *B. amyloliquefaciens* as well as *B. niger*, *B. brevis*, *B. megaterium* ), *Streptomyces*, *Lactobacillus* (especially *L. lactis*, *L. casei* ) and *Streptococcus* (especially *S. thermophilus*, *S. cremoris*, *S. agalactiae* ) come into question. Especially asporogenic strains of *Bacilli* or *Lactobacilli* are suitable.

Figure 3 shows the plasmid pTetLuc1 (SEQ ID NO: 1), in which the production of firefly luciferase (SEQ ID NO: 2) of *Photinus pyralis* (The gene encoding firefly luciferase was originally cloned and sequenced in the middle of the 1980's by DeWet, J. et al. (1987) Mol. Cell. Biol. 7, 725-737) can be switched on by the addition of a chemical agent belonging to the tetracycline family of antimicrobial agents in a cloned *E. coli* bacterium. SEQ ID NO: 1 shows the nucleotide sequence of this plasmid. By incubating *E. coli* cells containing this plasmid (or any other microbial strain where to similar regulation/reporter gene system is incorporated containing the necessary secondary regulatory sequences in the constructs such as correct ribosome binding region, transcriptional termination etc.) in the presence of

very small amounts of tetracyclines one is able to obtain light production by the addition of D-luciferin, which is the substrate of firefly luciferase. The intensity of light emission is proportional to the concentration of tetracycline used.

- 5 Figures 4a and 4b shows the effect of altogether seven different tetracyclines on the production of light as a function of concentration of each tetracycline. As controls different non-tetracycline antibiotics were included in this study to show that the sensor strain is specific for the tetracyclines. The luminescence was emitted from *E. coli* containing the plasmid pTetLux1 (SEQ ID NO: 3). The detection was made
- 10 after an incubation of 90 min. All tetracyclines tested behaved in a very similar manner and induction efficiencies were at the same antibiotic concentration area. This makes this sensor even more attractive for analytical use for the determination of the tetracycline group of antibiotics.
- 15 It should be noted that the accumulation of various tetracyclines into microbial cells is very strongly affected by the extracellular concentration of  $Mg^{2+}$  ions. Figure 5 shows the effect of increasing concentrations of  $Mg^{2+}$  ions on the behavior of *E. coli* cells containing the plasmid pTetLux1 (SEQ ID NO: 3). As can be seen the tetracycline response curve is shifted to the right as a function of increasing
- 20 concentrations of added  $Mg^{2+}$  ions. Thus by increasing the  $Mg^{2+}$  ion concentration one is able to decrease the sensitivity of the tetracycline sensor described in this invention. This fact is of great importance in cases where one does not need a high sensitivity of the measurement and where the approximate concentration of the ion is roughly constant and known such as in milk, serum and plasma.
- 25 The sensitivity can be increased by removing magnesium ions from the assay mixture e.g. by adding a chelating agent forming a complex with magnesium.

Figure 6 shows the possibilities to change the assay window for tetracyclines by adjusting the magnesium ion concentration and by combined adjustment of the magnesium ion concentration and pH.

- 5 The sensitivity of the assay can be increased by the use of cells which are especially antibiotic sensitive mutant strains. Hundreds of specific mutations for bacteria are known with which it is possible to study the activity of specific reactions. For instance trace amounts of antibiotics cause visible changes in the metabolism or in the cell membranes of antibiotic sensitive bacterial mutants. Mutations in cell wall
- 10 structural components or biosynthetic enzymes as well as in transport and efflux proteins such as porins might have an effect on the behavior of each sensor. Using these kinds of mutations one is able to develop tests measuring residual antibiotics from biological material very sensitively. It is also rather simple to transfer new characteristics into bacterial cells by genetic engineering techniques. This
- 15 phenomenon broadens the applicability of these organisms in tests utilizing whole cell sensor.

Measurement of light emission can be done by using X-ray or polaroid film, using a liquid scintillation counter, a CCD-camera or a luminometer. The CCD-camera is an instrument which is capable of detecting very low levels of light. In the applications of this invention such kind of a device could be used for the detection of tetracycline residues in food material such as vegetables or meat. The detection of light emission could be directly monitored from the surface of the food material sprayed with engineered luminescent bacteria. Either chemiluminescent (such as peroxidase -

25 luminol) or bioluminescent (such as luciferase - luciferin) reactions can be utilized. The luminometric method is performed with the aid of genes encoding either bacterial or beetle luciferases such as those described in the Figures 2 and 4. Several luminescent bacterial species such as *V. harveyi*, *V. fischeri*, *P. leiognathi*,

004240 29552560

SUB  
C2 20

- P. phosphoreum*, *Xenorhabdus luminescens* etc. exist. Luminescent beetles are for example *Luciola mingrelica*, *Photinus pyralis*, *Pyrophorus plagiophthalmus*, *Lampyris noctiluca*, *Pholas dactylus*, etc. Also several eukaryotic species in the sea which luminesce, such as marine ostracod *Vargula hilgendorfii*, jellyfish *Aequorea victoria*, batrachoidid fish *Porichtys notatus*, pempherid fish *Parapriacanthus ransonneti* etc. exist. Fluorescent reporter proteins such as green fluorescent protein (GFP) or any of its variants could be used in the methods described in this invention (Li, X. et al. (1997) J. Biol. Chem. 272, 28545-28549).
- 10 In this invention high detection sensitivity of the luminescent enzyme labels inside a living cell associated with tetracycline-specific induction of label synthesis is based on the use of optimal concentration of all the reactants inside the cell including the necessary cofactors and accessory enzymes. All luciferase genes from these organisms would presumably work in a similar manner as the two examples shown
- 15 in this invention. These systems together with enhancers and modulators (wavelength, emission kinetics etc.) of light emission has been described in more detail in Campbell, A. "Chemiluminescence; principles and applications in biology and medicine", Weinheim; Deerfield Beach, Fl.; VCH; Chichester: Horwood, 1988.
- 20 Peroxidases or oxidases can be used together with compounds such as luminol or acridines (for instance lucigenin) to yield luminescent signals suitable for a detection system described here. Enzymatically generated chemiluminescence offers great sensitivity and rapid detection, too, in assays described in this invention. Thermally stable dioxetanes (such as AMPPD and Lumigen PPD) can be
- 25 enzymatically (such as alkaline phosphatase or  $\beta$ -galactosidase) triggered to produce chemiluminescence (Schaap, A.P. et al. (1989) Clin. Chem. 35, 1863-1864). The only difference to the luciferase enzymes would be that these enzymes are capable

of cleaving a man-made substrate which gives light emission (chemiluminescence) and the luciferases cleave natural substrates to produce light (bioluminescence).

- Tetracycline-controlled expression systems are developed to express heterologous proteins in procaryotic and eucaryotic cells for the purpose of production under a tight control of tet-regulatory system ( Skerra, A. (1994) *Gene* 151, 131-135; Gossen, M. and Bujard, H. (1995) *US Patent* 5,464,758 ; Lutz, R. and Bujard, H. (1997) *Nucleic Acids Res.* 25, 1203-1210).
- 10 A method to study various tetracyclines and their mode of action was developed by Chopra et al. (Chopra, I. et al. (1990) *Antimicrob. Agents Chemother.* 34, 111-116) The assay system developed in this study was based on expression of  $\beta$ -galactosidase gene inserted under the control of tetA-gene. The method resulted in less sensitive detection of tetracyclines compared to the invention described here.
- 15 However in order to obtain maximum sensitivities Chopra et al. showed that it was necessary to add cyclic AMP (cAMP) to the medium which is an extremely expensive molecule to be used in routine applications. Furthermore, the method described by Chopra et al. contains a cell disruption stage by sonication in order to assay for the reporter gene activity,  $\beta$ -galactosidase, which step is not practical.
- 20 Instead, the method described in this invention does not contain any cell disruption. The activity of luciferase can be measured directly from living cells in real-time and in the case of pTetLux1 (SEQ ID NO: 3) there is no need of addition of any substrates. Therefore, promoter activation due to the presense/absense of tetracycline can be monitored continuously.

25

## EXPERIMENTS

As cloning hosts and in antibiotic residue measurements various *E. coli* MC1061 (cI+, *araD*139,  $\Delta$ (*ara-leu*)7696, *lacX*74, *galU*, *galK*, *hsr*, *hsm*, *strA*) (Casadaban,

M.J. and Cohen, S.N. (1980) J. Mol. Biol. 138, 179-207), BW322 (CGSC, *rfa210::Tn10*, *thi-1*, *relA1*, *spoT1*, *pyrE*) and K-12 (M72 Sm<sup>R</sup> *lacZm-ΔbiouvB*, *trpEA2*, Nam7Nam53cI857 HI) (Remaut, E. et al. (1981) Gene 15, 81-93) can be used. Especially the strain LH530 (Hirvas, L. et al. (1997) Microbiology 143, 73-81) which has a decreased rate of lipid A biosynthesis. It has proven to be hypersusceptible to many different antibiotics.

Cells were grown on appropriate minimal agar-plates and were kept maximally one month at +4 °C after which new plates were streaked. The strains were kept also in 15% glycerol at -70 °C, where from growth was started through minimal plates. The cells were first cultivated in 5 ml of 2xTY medium (16 g Bacto tryptone, 8 g Yeast extract, 8 g NaCl, H<sub>2</sub>O ad 1 l, pH 7.4, with appropriate antibiotic) 10 h at 30 °C in a shaker after which the cultivation was transferred to a bigger volume for 10 h with same medium.

#### Construction of tetracycline-responsive sensor plasmids:

To construct a recombinant DNA vector carrying luciferase genes under the control of a tetracycline responsive elements two new vectors were created. In the first one modified firefly luciferase gene (SEQ ID NO: 1) from vector pBLuc\* (Bonin, A.L. et al. (1994) Gene 141, 75-77) was excised by using restriction enzymes *XbaI* and *HinDIII* and the 1.7 kb fragment was isolated from LGT-agarose gel and purified using Qiagen gel extraction kit. This DNA-fragment containing the entire *Photinus pyralis* luciferase gene (SEQ ID NO: 1) was ligated using T4-DNA-ligase enzyme to vector pASK75 (Skerra, A. (1994) Gene 151, 131-135) which was previously restricted with the same restriction enzymes *XbaI* and *HinDIII* and calf intestinal phosphatase treated to remove the protruding phosphate groups in order to prevent self-ligation. The resulting ligation mixture was incubated 3 hours at room temperature after which one µl of the mixture was electroporated according to

Dower *et al.* (Dower, W.J. et al. (1988) Nucleic Acids Res. 16, 6126-6144) into electrocompetent *E. coli* MC1061 cells. A plasmid was extracted from one of the colonies obtained and checked for the estimated structure by appropriate restriction enzyme digestions and agarose gel electrophoretic techniques. The plasmid obtained  
5 was named as pTetLuc1 (SEQ ID NO: 1).

The plasmid containing the luxCDABE genes (SEQ ID NO: 3) of *Photobacterium luminescens* under the control of tetracycline responsive element was created as follows: Plasmid pASK75 was cut with restriction enzyme *EcoRI* and CIP-treated.  
10 The linearized plasmid was separated on a LGT-agarose gel electrophoresis and the agarose was removed by using the Qiagen kit. The lux operon was excised with *EcoRI* from plasmid pCGLS-11 (Frackman, S. et al. (1990) J. Bacteriol. 172, 5767-5773), gel purified as above and ligated to pASK75 by using T4-DNA-ligase at 16 °C overnight. The ligation mixture was electroporated into *E. coli* MC1061 cells as  
15 described above and correct transformants were screened for their ability to produce light (as measured with a BioOrbit 1250 manual luminometer) which production was increased in the presence of 1 µg/ml of tetracycline-HCl. The plasmid was further verified by restriction enzyme digestions and the correct structure was named as pTetLux1 (SEQ ID NO: 3). All the DNA-manipulations were performed  
20 according to Sambrook *et al.*, "Molecular Cloning: A laboratory Manual, Cold Spring Harbor Laboratory Press: Cold Spring Harbor, NY, 1989.

The vector pASK75 was utilized in the construction of tet-sensor plasmids shown in this invention. The vector pASK75 was originally developed for protein production  
25 and purification purposes. It contains a signal sequence for secretion of the recombinant protein into the periplasmic space of *E. coli*. Also a C-terminal fusion between a purification tail, the Strept-tag, was incorporated into the vector to facilitate purification of recombinant protein using streptavidin affinity agarose gel

09529967-042400

chromatography. The element controlling recombinant gene expression in the vector is tetA promoter/operator system that allows efficient regulation of the expression, which in Skerra's paper was described for the production and one-step purification of a murine single-chain antibody fragment. The tetA promoter/operator (SEQ ID NO: 9) is controlled by tetR-repressor (SEQ ID NO: 9) which is produced by the corresponding gene (SEQ ID NO: 9). Some of the above mentioned elements were eliminated from the present plasmids due to unnecessary features with respect to this invention.

10 **Transfer of the tetracycline sensor vectors to the antibiotic sensitive *E. coli* strain:**

Either pTetLux1 (SEQ ID NO: 3) or pTetLuc1 (SEQ ID NO: 1) was transformed into *E. coli* LH530 cells by electroporation as described above. The transformed cells were restreaked on agar plates and kept maximally for 2 weeks at +4 °C after which a new plate was streaked.

**Use of the manipulated *E. coli* in tetracycline determination methods:**

Example 1

Freeze-dried *E. coli* K-12/pTetLux1 were reconstituted with 1.0 ml of L-broth and bacteria were diluted 1:10 with 25 mM MES buffer in M9 minimal medium, pH 6.0. 190 µl bacterial suspension was added to microtiter plate wells containing 10 µl of tetracycline dilutions. The plate was incubated 90 minutes at 37 °C after which the plate was measured with Labsystems Luminoskan luminometer. As seen from Figure 7 the sensitivity of the assay of tetracycline is very high and comparable to that of fresh cells.



### Example 2

Two different types of sensor DNA vector construct were compared. Strains *E. coli* K-12/pTetLux1 and *E. coli* K-12/pTetLuc1 were cultivated in L-broth media until optical density measured at 600 nm (OD600) was 1.5. The cells were diluted 1 to 50 with 25 mM MES-buffer in M9 minimal medium, pH 6.0 (Sambrook *et al.*, 1989, Cold Spring Harbor Laboratory, Cold Spring Harbor) and 190 µl was added to microtitration plate wells and 10 µl of sample dilution of tetracycline was added. After a 60 min incubation at 37 °C the light emission was measured using a Labsystems Luminoskan luminometer. Figure 8 shows the bioluminescence dose response curve as a function of tetracycline added. As seen from the figure both systems (bacterial and insect luciferase) give roughly equal sensitivity of tetracycline detection.

One is able to use different luciferases instead of bacterial luciferase (SEQ ID NO: 4-8) from *P. luminescens* without losing sensitivity or other performance of the test. Figure 8 shows an analogous measurement to the one in Figure 4b. In the plasmid used in this test (pTetLuc1) the bacterial luciferase was compensated with firefly luciferase (SEQ ID NO: 2) as described in Figure 3. The test was done essentially as with bacterial luciferase except that after the cells had been incubated with or without tetracycline 10 minutes at 37 °C the cells were measured for light production after 15 minutes incubation time at 37 °C by adding 100 µl of solution containing 1 mM D-luciferin, in 0.1 M Na-citrate buffer, pH 5.0. Thereafter light production was measured using a manual luminometer 1250 (LKB-Wallac, Turku, Finland). As can be seen from Figure 8 sensitivity of the method to detect tetracycline hydrochloride is extremely high and comparable to the detection made with bacterial luciferase.

### Example 3

A lipemic pig serum was spiked at different concentrations of tetracycline, chlorotetracycline and oxytetracycline. Fresh *E. coli* K-12/pTetLux1 were diluted 1:50 with 25 mM MES buffer in M9 minimal medium, pH 6.0. 100 µl bacterial suspension was added to microtiter plate wells containing 100 µl of pig serum spiked with different tetracyclines. The plate was incubated 90 minutes at 37 °C after which the plate was measured with Labsystems Luminoskan luminometer. As seen from Figure 9 the sensitivity of the assay of different tetracyclines in pig serum matrix is very high.

10

### Example 4

Tetracyclines will form chelate complexes with  $\text{Ca}^{2+}$  and  $\text{Mg}^{2+}$  in samples (e.g. milk), and loose their antimicrobial and induction activity in our assay system. Tetracyclines can be displaced from cation chelates by using strong chelating agents such as EDTA. Figure 10 shows the determination of tetracycline from a milk sample, which is spiked with different concentrations of tetracycline. Different amounts of EDTA were added to milk samples and this kind of displacement of cation-tetracycline complex clearly improved the sensitivity of the assay. In the assay we used freeze-dried *E. coli* K12/pTetLux1 that were reconstituted with L-broth 10 minutes in room temperature before the assay.

15

20

### Example 5

Figure 11 shows the kinetics of bacterial bioluminescence after exposure of *E. coli* K-12/pTetLux1 to different dilutions of tetracycline antibiotics. The specific induction of tetracycline is very fast and specific light emission is seen already at the 10 minutes measuring point in the assay.

25

It will be appreciated that the methods of the present invention can be incorporated in the form of a variety of embodiments, only a few of which are disclosed herein. It will be apparent for the specialist in the field that other embodiments exist and do not depart from the spirit of the invention. Thus, the described embodiments are

5 illustrative and should not be construed as restrictive.

004240" 29662560

## SEQUENCE LISTING

## (1) GENERAL INFORMATION:

## (i) APPLICANT:

(A) NAME: KORPELA, Matti  
(B) STREET: Maijamaentie 13  
(C) CITY: Naantali  
(E) COUNTRY: Finland  
(F) POSTAL CODE (ZIP): FIN-21100

(A) NAME: KARP, Matti  
(B) STREET: Kampakatu 1  
(C) CITY: Kaarina  
(E) COUNTRY: Finland  
(F) POSTAL CODE (ZIP): FIN-20660

(A) NAME: KURITTU, Jussi  
(B) STREET: Puutarhakatu 16 A 20  
(C) CITY: Turku  
(E) COUNTRY: Finland  
(F) POSTAL CODE (ZIP): FIN-20100

(ii) TITLE OF INVENTION: A NEW ASSAY METHOD

(iii) NUMBER OF SEQUENCES: 11

## (iv) COMPUTER READABLE FORM:

(A) MEDIUM TYPE: Floppy disk  
(B) COMPUTER: IBM PC compatible  
(C) OPERATING SYSTEM: PC-DOS/MS-DOS  
(D) SOFTWARE: PatentIn Release #1.0, Version #1.30 (EPO)

## (vi) PRIOR APPLICATION DATA:

(A) APPLICATION NUMBER: FI 974235  
(B) FILING DATE: 14-NOV-1997

## (2) INFORMATION FOR SEQ ID NO: 1:

## (i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 4846 base pairs  
(B) TYPE: nucleic acid  
(C) STRANDEDNESS: double  
(D) TOPOLOGY: circular

(ii) MOLECULE TYPE: DNA (genomic)

(iii) HYPOTHETICAL: NO

(iv) ANTI-SENSE: NO

## (vi) ORIGINAL SOURCE:

(A) ORGANISM: Photinus pyralis

## (vii) IMMEDIATE SOURCE:

(B) CLONE: pTetLuc1

## (viii) POSITION IN GENOME:

(A) CHROMOSOME/SEGMENT: Plasmid

004240 29662560

## (ix) FEATURE:

- (A) NAME/KEY: misc\_feature  
(B) LOCATION:1..3098  
(D) OTHER INFORMATION:/standard\_name= "Vector pASK75"  
/note= "Part of plasmid originating from vector pASK75;  
feature description below, SEQ ID 9-11."  
/citation= ([2])

## (ix) FEATURE:

- (A) NAME/KEY: CDS  
(B) LOCATION:3119..4768  
(D) OTHER INFORMATION:/product= "Photinus pyralis  
luciferase"  
/citation= ([1])

## (x) PUBLICATION INFORMATION:

- (A) AUTHORS: Bonin,  
(B) TITLE: Photinus pyralis luciferase: vectors that  
contain a modified luc coding sequence allowing  
convenient transfer into other systems  
(C) JOURNAL: Gene  
(D) VOLUME: 141  
(F) PAGES: 75-77  
(G) DATE: 1994  
(K) RELEVANT RESIDUES IN SEQ ID NO: 1: FROM 3099 TO 4772

## (x) PUBLICATION INFORMATION:

- (A) AUTHORS: Skerra, A  
(B) TITLE: Use of the tetracycline promoter for the  
tightly regulated production of a murine antibody  
fragment in Escherichia coli  
(C) JOURNAL: Gene  
(D) VOLUME: 151  
(E) ISSUE: 1-2  
(F) PAGES: 131-135  
(G) DATE: 30-DEC-1994  
(K) RELEVANT RESIDUES IN SEQ ID NO: 1: FROM 1 TO 3098

## (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 1:

AGCTTGACCT GTGAAGTGAA AAATGGCGCA CATTGTGCGA CATTTTTTTTT GTCTGCCGTT 60  
TACCGCTACT GCGTCACGGA TCTCCACGCG CCCTGTAGCG GCGCATTAAAG CGCGGCGGGT 120  
GTGGTGGTTA CGCGCAGCGT GACCGCTACA CTTGCCAGCG CCCTAGCGCC CGCTCCTTTC 180  
GCTTCTTCC CTTCTTTTCT CGCCACGTTT GCCGGCTTTC CCCGTCAAGC TCTAAATCGG 240  
GGGCTCCCTT TAGGGTTCCG ATTTAGTGCT TTACGGCACC TCGACCCCAA AAAACTTGAT 300  
TAGGGTGATG GTTACGCTAG TGGGCCATCG CCCTGATAGA CGGTTTTTCG CCCTTTGACG 360  
TTGGAGTCCA CGTTCTTTAA TAGTGGACTC TTGTTCCAAA CTGGAACAAC ACTCAACCCT 420  
ATCTGGGTCT ATTCTTTTGA TTTATAAGGG ATTTTGCCGA TTTCGGCCTA TTGGTTAAAA 480  
AATGAGCTGA TTTAACAAAA ATTTAACGCG AATTTTAACA AAATATTAAC GCTTACAATT 540  
TCAGGTGGCA CTTTTCGGGG AAATGTGCGC GGAACCCCTA TTTGTTTATT TTTCTAAATA 600  
CATTCAAATA TGTATCCGCT CATGAGACAA TAACCCTGAT AAATGCTTCA ATAATATTGA 660  
AAAAGGAAGA GTATGAGTAT TCAACATTTT CGTGTCGCCC TTATTCCCTT TTTTGC GGCA 720  
TTTTGCCTTC CTGTTTTTGC TCACCCAGAA ACGCTGGTGA AAGTAAAGA TGCTGAAGAT 780  
CAGTTGGGTG CACGAGTGGG TTACATCGAA CTGGATCTCA ACAGCGGTAA GATCCTTGAG 840

004240 29662560

AGTTTTTCGCC CCGAAGAACG TTTTCCAATG ATGAGCACTT TTAAAGTTCT GCTATGTGGC 900  
GCGGTATTAT CCCGTATTGA CGCCGGGCAA GAGCAACTCG GTCGCCGCAT AACTATTCT 960  
CAGAACTGACT TGGTTGAGTA CTCACCAGTC ACAGAAAAGC ATCTTACGGA TGGCATGACA 1020  
GTAAGAGAAT TATGCAGTGC TGCCATAACC ATGAGTGATA AACTGCGGC CAACTTACTT 1080  
CTGACAACGA TCGGAGGACC GAAGGAGCTA ACCGCTTTTT TGCACAACAT GGGGGATCAT 1140  
GTAACTCGCC TTGATCGTTG GGAACCGGAG CTGAATGAAG CCATACCAA CGACGAGCGT 1200  
GACACCACGA TGCCTGTAGC AATGGCAACA ACGTTGCGCA AACTATTAA TGGCGAACTA 1260  
CTTACTCTAG CTTCCCGGCA ACAATTGATA GACTGGATGG AGGCGGATAA AGTTGCAGGA 1320  
CCACTTCTGC GCTCGGCCCT TCCGGCTGGC TGGTTTATTG CTGATAAATC TGGAGCCGGT 1380  
GAGCGTGGCT CTCGCGGTAT CATTGCAGCA CTGGGGCCAG ATGGTAAGCC CTCCCGTATC 1440  
GTAGTTATCT ACACGACGGG GAGTCAGGCA ACTATGGATG AACGAAATAG ACAGATCGCT 1500  
GAGATAGGTG CCTCACTGAT TAAGCATTGG TAGGAATTAA TGATGTCTCG TTTAGATAAA 1560  
AGTAAAGTGA TTAACAGCGC ATTAGAGCTG CTTAATGAGG TCGGAATCGA AGGTTTAACA 1620  
ACCCGTAAAC TCGCCAGAA GCTAGGTGTA GAGCAGCCTA CATTGTATTG GCATGTAAAA 1680  
AATAAGCGGG CTTTGCTCGA CGCCTTAGCC ATTGASATGT TAGATAGGCA CCATACTCAC 1740  
TTTTGCCCCT TAGAAGGGGA AAGCTGGCAA GATTTTTTAC GTAATAACGC TAAAAGTTTT 1800  
AGATGTGCTT TACTAAGTCA TCGCGATGGA GCAAAAGTAC ATTTAGGTAC ACGGCCTACA 1860  
GAAAAACAGT ATGAACTCT CGAAAAATCAA TTAGCCTTTT TATGCCAACA AGGTTTTTCA 1920  
CTAGAGAATG CATTATATGC ACTCAGCGCA GTGGGGCATT TTACTTTAGG TTGCGTATTG 1980  
GAAGATCAAG AGCATCAAGT CGTTAAAGAA GAAAGGGAAA CACCTACTAC TGATAGTATG 2040  
CCGCCATTAT TACGACAAGC TATCGAATTA TTTGATCACC AAGGTGCAGA GCCAGCCTTC 2100  
TTATTCGGCC TTGAATTGAT CATATGCGGA TTAGAAAAAC AACTTAAATG TGAAAGTGGG 2160  
TCTTAAAGC AGCATAACCT TTTCCGTGA TGGTAACTTC ACTAGTTTAA AAGGATCTAG 2220  
GTGAAGATCC TTTTGTGATA TCTCATGACC AAAATCCCTT AACGTGAGTT TTCGTTCCAC 2280  
TGAGCGTCAG ACCCGTAGA AAAGATCAA GGATCTTCTT GAGATCCTTT TTTCTGCGC 2340  
GTAATCTGCT GCTTGCAAAC AAAAAACCA CCGTACCAG CGGTGGTTTG TTTGCCGGAT 2400  
CAAGAGCTAC CAACTCTTTT TCCGAAGGTA ACTGGCTTCA GCAGAGCGCA GATACCAAAT 2460  
ACTGTCCTTC TAGTGTAGCC GTAGTTAGGC CACCACTTCA AGAACTCTGT AGCACCGCCT 2520  
ACATACCTCG CTCTGCTAAT CCTGTTACCA GTGGCTGCTG CCAGTGGCGA TAAGTCGTGT 2580  
CTTACCGGGT TGGACTCAAG ACGATAGTTA CCGGATAAGG CGCAGCGGTC GGGCTGAACG 2640  
GGGGTTTCGT GCACACAGCC CAGCTTGGAG CGAACGACCT ACACCGAACT GAGATACCTA 2700  
CAGCGTGAGC TATGAGAAAG CGCCACGCTT CCCGAAGGGA GAAAGGCGGA CAGGTATCCG 2760  
GTAAGCGGCA GGGTCGGAAC AGGAGAGCGC ACGAGGGAGC TTCCAGGGGG AAACGCCTGG 2820  
TATCTTTATA GTCCTGTCGG GTTTCGCCAC CTCTGACTTG AGCGTCGATT TTTGTGATGC 2880

092997 042400

TCGTCAGGGG GCGGAGCCT ATGGA AAAAC GCCAGCAACG CGGCCTTTT ACGGTTCCTG 2940  
 GCCTTTTGCT GGCCTTTTGC TCACATGACC CGACACCATC GAATGGCCAG ATGATTAATT 3000  
 CCTAATTTTT GTTGACACTC TATCATTGAT AGAGTTATTT TACCAC TCCC TATCAGTGAT 3060  
 AGAGAAAAGT GAAATGAATA GTTCGACAAA AATCTAGAAC TAGTGGATCC CCCGTACC 3118  
 ATG GAA GAC GCC AAA AAC ATA AAG AAA GGC CCG GCG CCA TTC TAT CCG 3166  
 Met Glu Asp Ala Lys Asn Ile Lys Lys Gly Pro Ala Pro Phe Tyr Pro  
 1 5 10 15  
 CTA GAG GAT GGA ACC GCT GGA GAG CAA CTG CAT AAG GCT ATG AAG AGA 3214  
 Leu Glu Asp Gly Thr Ala Gly Glu Gln Leu His Lys Ala Met Lys Arg  
 20 25 30  
 TAC GCC CTG GTT CCT GGA ACA ATT GCT TTT ACA GAT GCA CAT ATC GAG 3262  
 Tyr Ala Leu Val Pro Gly Thr Ile Ala Phe Thr Asp Ala His Ile Glu  
 35 40 45  
 GTG AAC ATC ACG TAC GCG GAA TAC TTC GAA ATG TCC GTT CGG TTG GCA 3310  
 Val Asn Ile Thr Tyr Ala Glu Tyr Phe Glu Met Ser Val Arg Leu Ala  
 50 55 60  
 GAA GCT ATG AAA CGA TAT GGG CTG AAT ACA AAT CAC AGA ATC GTC GTA 3358  
 Glu Ala Met Lys Arg Tyr Gly Leu Asn Thr Asn His Arg Ile Val Val  
 65 70 75 80  
 TGC AGT GAA AAC TCT CTT CAA TTC TTT ATG CCG GTG TTG GGC GCG TTA 3406  
 Cys Ser Glu Asn Ser Leu Gln Phe Phe Met Pro Val Leu Gly Ala Leu  
 85 90 95  
 TTT ATC GGA GTT GCA GTT GCG CCC GCG AAC GAC ATT TAT AAT GAA CGT 3454  
 Phe Ile Gly Val Ala Val Ala Pro Ala Asn Asp Ile Tyr Asn Glu Arg  
 100 105 110  
 GAA TTG CTC AAC AGT ATG AAC ATT TCG CAG CCT ACC GTA GTG TTT GTT 3502  
 Glu Leu Leu Asn Ser Met Asn Ile Ser Gln Pro Thr Val Val Phe Val  
 115 120 125  
 TCC AAA AAG GGG TTG CAA AAA ATT TTG AAC GTG CAA AAA AAA TTA CCA 3550  
 Ser Lys Lys Gly Leu Gln Lys Ile Leu Asn Val Gln Lys Lys Leu Pro  
 130 135 140  
 ATA ATC CAG AAA ATT ATT ATC ATG GAT TCT AAA ACG GAT TAC CAG GGA 3598  
 Ile Ile Gln Lys Ile Ile Ile Met Asp Ser Lys Thr Asp Tyr Gln Gly  
 145 150 155 160  
 TTT CAG TCG ATG TAC ACG TTC GTC ACA TCT CAT CTA CCT CCC GGT TTT 3646  
 Phe Gln Ser Met Tyr Thr Phe Val Thr Ser His Leu Pro Pro Gly Phe  
 165 170 175  
 AAT GAA TAC GAT TTT GTA CCA GAG TCC TTT GAT CGT GAC AAA ACA ATT 3694  
 Asn Glu Tyr Asp Phe Val Pro Glu Ser Phe Asp Arg Asp Lys Thr Ile  
 180 185 190  
 GCA CTG ATA ATG AAC TCC TCT GGA TCT ACT GGG TTA CCT AAG GGT GTG 3742  
 Ala Leu Ile Met Asn Ser Ser Gly Ser Thr Gly Leu Pro Lys Gly Val  
 195 200 205  
 GCC CTT CCG CAT AGA ACT GCC TGC GTC AGA TTC TCG CAT GCC AGA GAT 3790  
 Ala Leu Pro His Arg Thr Ala Cys Val Arg Phe Ser His Ala Arg Asp  
 210 215 220  
 CCT ATT TTT GGC AAT CAA ATC ATT CCG GAT ACT GCG ATT TTA AGT GTT 3838  
 Pro Ile Phe Gly Asn Gln Ile Ile Pro Asp Thr Ala Ile Leu Ser Val  
 225 230 235 240

001210 29662560

GTT CCA TTC CAT CAC GGT TTT GGA ATG TTT ACT ACA CTC GGA TAT TTG 3886  
 Val Pro Phe His His Gly Phe Gly Met Phe Thr Thr Leu Gly Tyr Leu 255  
 245 250

ATA TGT GGA TTT CGA GTC GTC TTA ATG TAT AGA TTT GAA GAA GAG CTG 3934  
 Ile Cys Gly Phe Arg Val Val Leu Met Tyr Arg Phe Glu Glu Glu Leu 270  
 260 265

TTT TTA CGA TCC CTT CAG GAT TAC AAA ATT CAA AGT GCG TTG CTA GTA 3982  
 Phe Leu Arg Ser Leu Gln Asp Tyr Lys Ile Gln Ser Ala Leu Leu Val 285  
 275 280

CCA ACC CTA TTT TCA TTC TTC GCC AAA AGC ACT CTG ATT GAC AAA TAC 4030  
 Pro Thr Leu Phe Ser Phe Phe Ala Lys Ser Thr Leu Ile Asp Lys Tyr 300  
 290 295

GAT TTA TCT AAT TTA CAC GAA ATT GCT TCT GGG GGC GCA CCT CTT TCG 4078  
 Asp Leu Ser Asn Leu His Glu Ile Ala Ser Gly Gly Ala Pro Leu Ser 320  
 305 310 315

AAA GAA GTC GGG GAA GCG GTT GCA AAA CGC TTC CAT CTT CCA GGG ATA 4126  
 Lys Glu Val Gly Glu Ala Val Ala Lys Arg Phe His Leu Pro Gly Ile 335  
 325 330

CGA CAA GGA TAT GGG CTC ACT GAG ACT ACA TCA GCT ATT CTG ATT ACA 4174  
 Arg Gln Gly Tyr Gly Leu Thr Glu Thr Thr Ser Ala Ile Leu Ile Thr 350  
 340 345

CCC GAG GGG GAT GAT AAA CCG GGC GCG GTC GGT AAA GTT GTT CCA TTT 4222  
 Pro Glu Gly Asp Asp Lys Pro Gly Ala Val Gly Lys Val Val Pro Phe 365  
 355 360

TTT GAA GCG AAG GTT GTG GAT CTG GAT ACC GGG AAA ACG CTG GGC GTT 4270  
 Phe Glu Ala Lys Val Val Asp Leu Asp Thr Gly Lys Thr Leu Gly Val 380  
 370 375

AAT CAG AGA GGC GAA TTA TGT GTC AGA GGA CCT ATG ATT ATG TCC GGT 4318  
 Asn Gln Arg Gly Glu Leu Cys Val Arg Gly Pro Met Ile Met Ser Gly 400  
 385 390 395

TAT GTA AAC AAT CCG GAA GCG ACC AAC GCC TTG ATT GAC AAG GAT GGA 4366  
 Tyr Val Asn Asn Pro Glu Ala Thr Asn Ala Leu Ile Asp Lys Asp Gly 415  
 405 410

TGG CTA CAT TCT GGA GAC ATA GCT TAC TGG GAC GAA GAC GAA CAC TTC 4414  
 Trp Leu His Ser Gly Asp Ile Ala Tyr Trp Asp Glu Asp Glu His Phe 430  
 420 425

TTC ATA GTT GAC CGC TTG AAG TCT TTA ATT AAA TAC AAA GGA TAC CAG 4462  
 Phe Ile Val Asp Arg Leu Lys Ser Leu Ile Lys Tyr Lys Gly Tyr Gln 445  
 435 440

GTG GCC CCC GCT GAA TTG GAG TCG ATA TTG TTA CAA CAC CCC AAC ATC 4510  
 Val Ala Pro Ala Glu Leu Glu Ser Ile Leu Leu Gln His Pro Asn Ile 460  
 450 455

TTC GAC GCG GGC GTG GCA GGT CTT CCC GAC GAT GAC GCC GGT GAA CTT 4558  
 Phe Asp Ala Gly Val Ala Gly Leu Pro Asp Asp Asp Ala Gly Glu Leu 480  
 465 470 475

CCC GTC GCC GTT GTT GTT TTG GAG CAC GGA AAG ACG ATG ACG GAA AAA 4606  
 Pro Ala Ala Val Val Val Leu Glu His Gly Lys Thr Met Thr Glu Lys 495  
 485 490

GAG ATC GTG GAT TAC GTC GCC AGT CAA GTA ACA ACC GCC AAA AAG TTG 4654  
 Glu Ile Val Asp Tyr Val Ala Ser Gln Val Thr Thr Ala Lys Lys Leu 510  
 500 505

004240" 29662560



CGC GGA GGA GTT GTG TTT GTG GAC GAA GTA CCG AAA GGT CTT ACC GGA 4702  
 Arg Gly Gly Val Val Phe Val Asp Glu Val Pro Lys Gly Leu Thr Gly  
 515 520 525

AAA CTC GAC GCA AGA AAA ATC AGA GAG ATC CTC ATA AAG GCC AAG AAG 4750  
 Lys Leu Asp Ala Arg Lys Ile Arg Glu Ile Leu Ile Lys Ala Lys Lys  
 530 535 540

GGC GGA AAG TCC AAA TTG TAAAATGTAA CTGTATTTCAG CGATGACGAA 4798  
 Gly Gly Lys Ser Lys Leu  
 545 550

ATTCTTAGCT ATTGTAATAC TCTAGCGGGC TGCAGGAATT CGATATCA 4846

## (2) INFORMATION FOR SEQ ID NO: 2:

- (i) SEQUENCE CHARACTERISTICS:  
 (A) LENGTH: 550 amino acids  
 (B) TYPE: amino acid  
 (D) TOPOLOGY: linear

- (ii) MOLECULE TYPE: protein  
 (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 2:

Met Glu Asp Ala Lys Asn Ile Lys Lys Gly Pro Ala Pro Phe Tyr Pro  
 1 5 10 15

Leu Glu Asp Gly Thr Ala Gly Glu Gln Leu His Lys Ala Met Lys Arg  
 20 25 30

Tyr Ala Leu Val Pro Gly Thr Ile Ala Phe Thr Asp Ala His Ile Glu  
 35 40 45

Val Asn Ile Thr Tyr Ala Glu Tyr Phe Glu Met Ser Val Arg Leu Ala  
 50 55 60

Glu Ala Met Lys Arg Tyr Gly Leu Asn Thr Asn His Arg Ile Val Val  
 65 70 75 80

Cys Ser Glu Asn Ser Leu Gln Phe Phe Met Pro Val Leu Gly Ala Leu  
 85 90 95

Phe Ile Gly Val Ala Val Ala Pro Ala Asn Asp Ile Tyr Asn Glu Arg  
 100 105 110

Glu Leu Leu Asn Ser Met Asn Ile Ser Gln Pro Thr Val Val Phe Val  
 115 120 125

Ser Lys Lys Gly Leu Gln Lys Ile Leu Asn Val Gln Lys Lys Leu Pro  
 130 135 140

Ile Ile Gln Lys Ile Ile Ile Met Asp Ser Lys Thr Asp Tyr Gln Gly  
 145 150 155 160

Phe Gln Ser Met Tyr Thr Phe Val Thr Ser His Leu Pro Pro Gly Phe  
 165 170 175

Asn Glu Tyr Asp Phe Val Pro Glu Ser Phe Asp Arg Asp Lys Thr Ile  
 180 185 190

Ala Leu Ile Met Asn Ser Ser Gly Ser Thr Gly Leu Pro Lys Gly Val  
 195 200 205

Ala Leu Pro His Arg Thr Ala Cys Val Arg Phe Ser His Ala Arg Asp  
 210 215 220

004240 29562960

Pro Ile Phe Gly Asn Gln Ile Ile Pro Asp Thr Ala Ile Leu Ser Val  
225 230 235 240  
Val Pro Phe His His Gly Phe Gly Met Phe Thr Thr Leu Gly Tyr Leu  
245 250 255  
Ile Cys Gly Phe Arg Val Val Leu Met Tyr Arg Phe Glu Glu Glu Leu  
260 265 270  
Phe Leu Arg Ser Leu Gln Asp Tyr Lys Ile Gln Ser Ala Leu Leu Val  
275 280 285  
Pro Thr Leu Phe Ser Phe Phe Ala Lys Ser Thr Leu Ile Asp Lys Tyr  
290 295 300  
Asp Leu Ser Asn Leu His Glu Ile Ala Ser Gly Gly Ala Pro Leu Ser  
305 310 315 320  
Lys Glu Val Gly Glu Ala Val Ala Lys Arg Phe His Leu Pro Gly Ile  
325 330 335  
Arg Gln Gly Tyr Gly Leu Thr Glu Thr Thr Ser Ala Ile Leu Ile Thr  
340 345 350  
Pro Glu Gly Asp Asp Lys Pro Gly Ala Val Gly Lys Val Val Pro Phe  
355 360 365  
Phe Glu Ala Lys Val Val Asp Leu Asp Thr Gly Lys Thr Leu Gly Val  
370 375 380  
Asn Gln Arg Gly Glu Leu Cys Val Arg Gly Pro Met Ile Met Ser Gly  
385 390 395 400  
Tyr Val Asn Asn Pro Glu Ala Thr Asn Ala Leu Ile Asp Lys Asp Gly  
405 410 415  
Trp Leu His Ser Gly Asp Ile Ala Tyr Trp Asp Glu Asp Glu His Phe  
420 425 430  
Phe Ile Val Asp Arg Leu Lys Ser Leu Ile Lys Tyr Lys Gly Tyr Gln  
435 440 445  
Val Ala Pro Ala Glu Leu Glu Ser Ile Leu Leu Gln His Pro Asn Ile  
450 455 460  
Phe Asp Ala Gly Val Ala Gly Leu Pro Asp Asp Asp Ala Gly Glu Leu  
465 470 475 480  
Pro Ala Ala Val Val Val Leu Glu His Gly Lys Thr Met Thr Glu Lys  
485 490 495  
Glu Ile Val Asp Tyr Val Ala Ser Gln Val Thr Thr Ala Lys Lys Leu  
500 505 510  
Arg Gly Gly Val Val Phe Val Asp Glu Val Pro Lys Gly Leu Thr Gly  
515 520 525  
Lys Leu Asp Ala Arg Lys Ile Arg Glu Ile Leu Ile Lys Ala Lys Lys  
530 535 540  
Gly Gly Lys Ser Lys Leu  
545 550

004240" 29562560

## (2) INFORMATION FOR SEQ ID NO: 3:

- (i) SEQUENCE CHARACTERISTICS:  
    (A) LENGTH: 10220 base pairs  
    (B) TYPE: nucleic acid  
    (C) STRANDEDNESS: double  
    (D) TOPOLOGY: circular
- (ii) MOLECULE TYPE: DNA (genomic)
- (iii) HYPOTHETICAL: NO
- (iv) ANTI-SENSE: NO
- (vi) ORIGINAL SOURCE:  
    (A) ORGANISM: Photorhabdus luminescens
- (vii) IMMEDIATE SOURCE:  
    (B) CLONE: pTetLux1
- (ix) FEATURE:  
    (A) NAME/KEY: misc\_feature  
    (B) LOCATION: join(1..3190, 10140..10220)  
    (D) OTHER INFORMATION: /standard\_name= "vector pASK75"  
        /note= "Parts of plasmid originating from vector pASK75;  
            feature description below, SEQ ID NO: 9-11."  
        /citation= ([2])
- (ix) FEATURE:  
    (A) NAME/KEY: CDS  
    (B) LOCATION: 3634..5082  
    (D) OTHER INFORMATION: /product= "Lux C"  
        /citation= ([1])
- (ix) FEATURE:  
    (A) NAME/KEY: CDS  
    (B) LOCATION: 5097..6017  
    (D) OTHER INFORMATION: /product= "Lux D"  
        /citation= ([1])
- (ix) FEATURE:  
    (A) NAME/KEY: CDS  
    (B) LOCATION: 6069..7148  
    (D) OTHER INFORMATION: /product= "Lux A"  
        /citation= ([1])
- (ix) FEATURE:  
    (A) NAME/KEY: CDS  
    (B) LOCATION: 7166..8146  
    (D) OTHER INFORMATION: /product= "Lux B"  
        /citation= ([1])
- (ix) FEATURE:  
    (A) NAME/KEY: CDS  
    (B) LOCATION: 8256..9437  
    (D) OTHER INFORMATION: /product= "Lux E"  
        /citation= ([1])
- (x) PUBLICATION INFORMATION:  
    (A) AUTHORS: Frackman,  
    (B) TITLE: Cloning, organization and expression of the  
        bioluminescence genes of Xenorhabdus  
        luminescens  
    (C) JOURNAL: J. Bacteriol.  
    (D) VOLUME: 172  
    (F) PAGES: 5767-5773  
    (G) DATE: 1990  
    (K) RELEVANT RESIDUES IN SEQ ID NO: 3: FROM 3191 TO 10139

00529067.042400

## (x) PUBLICATION INFORMATION:

- (A) AUTHORS: Skerra, A  
(B) TITLE: Use of the tetracycline promoter for the tightly regulated production of a murine antibody fragment in Escherichia coli  
(C) JOURNAL: Gene  
(D) VOLUME: 151  
(E) ISSUE: 1-2  
(F) PAGES: 131-135  
(G) DATE: 30-DEC-1994  
(K) RELEVANT RESIDUES IN SEQ ID NO: 3: FROM 1 TO 3190

## (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 3:

AGCTTGACCT GTGAAGTGAA AAATGGCGCA CATTGTGCGA CATTTTTTTT GTCTGCCGTT 60  
TACCGCTACT GCGTCACGGA TCTCCACGCG CCCTGTAGCG GCGCATTAAAG CGCGGCGGGT 120  
GTGGTGGTTA CGCGCAGCGT GACCGCTACA CTTGCCAGCG CCCTAGCGCC CGCTCCTTTC 180  
GCTTTCTTCC CTTCCTTTCT CGCCACGTTT GCCGCTTTC CCCGTCAAGC TCTAAATCGG 240  
GGGCTCCCTT TAGGGTTCCG ATTTAGTGCT TTACGGCACC TCGACCCCAA AAAACTTGAT 300  
TAGGGTGATG GTTCACGTAG TGGGCCATCG CCCTGATAGA CGGTTTTTCG CCTTTTGACG 360  
TTGGAGTCCA CGTTCTTTAA TAGTGGACTC TTGTCCAAA CTGGAACAAC ACTCAACCCT 420  
ATCTCGGTCT ATTCTTTTGA TTTATAAGGG ATTTTGCCGA TTTCGGCCTA TTGGTTAAAA 480  
AATGAGCTGA TTTAACAAAA ATTTAAACCG AATTTTAACA AAATATTAAC GCTTACAATT 540  
TCAGGTGGCA CTTTTCGGGG AAATGTGGGC GGAACCCCTA TTTGTTTATT TTTCTAAATA 600  
CATTCAAATA TGTATCCGCT CATGAGACAA TAACCCTGAT AAATGCTTCA ATAATATTGA 660  
AAAAGGAAGA GTATGAGTAT TCAACATTTC CGTGTCGCCC TTATTCCTT TTTTGC GGCA 720  
TTTTGCCTTC CTGTTTTTGC TCACCCAGAA ACGCTGGTGA AAGTAAAAGA TGCTGAAGAT 780  
CAGTTGGGTG CACGAGTGGG TTACATCGAA CTGGATCTCA ACAGCGGTAA GATCCTTGAG 840  
AGTTTTTCGCC CCGAAGAACG TTTTCCAATG ATGAGCACTT TTAAAGTTCT GCTATGTGGC 900  
GCGGTATTAT CCCGTATTGA CGCCGGGCAA GAGCAACTCG GTCGCCGCAT ACACTATTCT 960  
CAGAATGACT TGGTTGAGTA CTCACCAGTC ACAGAAAAGC ATCTTACGGA TGGCATGACA 1020  
GTAAGAGAAT TATGCAGTGC TGCCATAACC ATGAGTGATA ACACTGCGGC CAACTTACTT 1080  
CTGACAACGA TCGGAGGACC GAAGGAGCTA ACCGCTTTTT TGCACAACAT GGGGGATCAT 1140  
GTAACTCGCC TTGATCGTTG GGAACCGGAG CTGAATGAAG CCATACCAAA CGACGAGCGT 1200  
GACACCAOGA TGCCTGTAGC AATGGCAACA ACGTTGCGCA AACTATTAAAC TGGCGAACTA 1260  
CTTACTCTAG CTTCCCGGCA ACAATTGATA GACTGGATGG AGGCGGATAA AGTTGCAGGA 1320  
CCACTCTGCT GCTCGGCCCT TCCGGCTGGC TGGTTTATTG CTGATAAATC TGGAGCCGGT 1380  
GAGCGTGGCT CTCGCGGTAT CATTGCAGCA CTGGGGCCAG ATGGTAAGCC CTCCCGTATC 1440  
GTAGTTATCT ACACGACGGG GAGTCAGGCA ACTATGGATG AACGAAATAG ACAGATCGCT 1500  
GAGATAGGTG CCTCACTGAT TAAGCATTGG TAGGAATTAA TGATGTCTCG TTTAGATAAA 1560

0952996 04240

AGTAAAGTGA TTAACAGCGC ATTAGAGCTG CTTAATGAGG TCGGAATCGA AGGTTTAAACA 1620  
ACCCGTAAAC TCGCCCAGAA GCTAGGTGTA GAGCAGCCTA CATTGTATTG GCATGTAAAA 1680  
AATAAGCGGG CTTTGCTCGA CGCCTTAGCC ATTGAGATGT TAGATAGGCA CCATACTCAC 1740  
TTTGGCCCTT TAGAAGGGGA AAGCTGGCAA GATTTTITAC GTAATAACGC TAAAAGTTT 1800  
AGATGTGCTT TACTAAGTCA TCGCGATGGA GCAAAAGTAC ATTTAGGTAC ACGGCCTACA 1860  
GAAAAACAGT ATGAAACTCT CGAAATCAA TTAGCCTTTT TATGCCAACA AGGTTTTTCA 1920  
CTAGAGAATG CATTATATGC ACTCAGCGCA GTGGGGCATT TTACTIONTAGG TTGCGTATTG 1980  
GAAGATCAAG AGCATCAAGT CGCTAAAGAA GAAAGGGAAA CACCTACTAC TGATAGTATG 2040  
CCGCCATTAT TACGACAAGC TATCGAATTA TTTGATCACC AAGGTGCAGA GCCAGCCTTC 2100  
TTATTCGGCC TTGAATTGAT CATATGCGGA TTAGAAAAAC AACTTAAATG TGAAAGTGGG 2160  
TCTTAAAGC AGCATAACCT TTTTCCGTGA TGTAACCTC ACTAGTTTAA AAGGATCTAG 2220  
GTGAAGATCC TTTTGTATA TCTCATGACC AAAATCCCTT AACGTGAGTT TTCGTTCCAC 2280  
TGAGCGTCAG ACCCCGTAGA AAAATCAAAA GGATCTTCTT GAGATCCTTT TTTCTGCGC 2340  
GTAATCTGCT GCTTGCAAAC AAAAABACCA CCGCTACCAG CGGTGGTTTG TTTGCCGGAT 2400  
CAAGAGCTAC CAACTCTTTT TCCGAGGTA ACTGGCTTCA GCAGAGCGCA GATACCAAAT 2460  
ACTGTCTTC TAGTGTAGCC GTAGTTAGGC CACCCTTCA AGAACTCTGT AGCACCAGCT 2520  
ACATACCTCG CTCTGCTAAT CCTGTTACCA GTGGCTGCTG CCAGTGGCGA TAAGTCGTGT 2580  
CTTACCGGGT TGGACTCAAG ACGATAGTTA CCGGATAAGG CGCAGCGGTC GGGCTGAACG 2640  
GGGGGTTTGT GCACACAGCC CAGCTTGGAG CGAACGACCT ACACCGAACT GAGATACCTA 2700  
CAGCGTGAGC TATGAGAAAAG CGCCACGCTT CCCGAAGGGA GAAAGGCGGA CAGGTATCCG 2760  
GTAAGCGGCA GGGTCGGAAC AGGAGAGCGC ACGAGGGAGC TTCCAGGGGG AAACGCCTGG 2820  
TATCTTTATA GTCTGTGCGG GTTTCGCCAC CTCTGACTTG AGCGTCGATT TTTGTGATGC 2880  
TCGTCAGGGG GCGCGAGCCT ATGGAAAAAC GCCAGCAACG CGGCCTTTTT ACGGTTCTCTG 2940  
GCCTTTTGTG GGCCTTTTGC TCACATGACC CGACACCATC GAATGGCCAG ATGATTAATT 3000  
CCTAATTTTT GTTGACACTC TATCATGAT AGAGTTATTT TACCACTCCC TATCAGTGAT 3060  
AGAGAAAAGT GAAATGAATA GTTCGACAAA AATCTAGATA ACGAGGGCAA AAAATGAAAA 3120  
AGACAGCTAT CGCGATTGCA GTGGCACTGG CTGGTTTCGC TACCGTAGCG CAGGCCTGAG 3180  
ACCAGAATTC TTCTTTAGAA ATCTGCCGGT AAAAAATTAGA TTGCTATTCA ATCTATTTCT 3240  
ATCGTATTT GTGAAATAAT ACTCAGGATA ATAATTTACA TAAATATTAT CACGCATTAG 3300  
AGAAGAGCAT GACTTTTTTA ATTTAACTT TTCATTAACA AATCTTGTTG ATATGAAAAT 3360  
TTTCCTTTGC TATTTTAAACA GATATTAAAA CGGGAATAGG CGTTATATTG ACGATCCATT 3420  
CAGTTAGATT AAAACCTTG AGCAGAAAAT TTATATTATT ATCATAATTA TGACGAAAGT 3480  
TACAGGCCAG GAACCACGTA GTCAGAATCT GATTTTCTAT ATATTTGTTA TTTACATCGT 3540  
CATAACACAA AAATATAAGA AGCAAGTGTT GGTACGACCA GTTCGCAAGA TAGTTAAACA 3600

004240 2962560

GCAACTTAAG TTGAAATTAC CCCCATTAAG TGG ATG GCA AAT ATG ACT AAA AAA 3654  
 Met Ala Asn Met Thr Lys Lys 555

ATT TCA TTC ATT ATT AAC GGC CAG GTT GAA ATC TTT CCC GAA AGT GAT 3702  
 Ile Ser Phe Ile Ile Asn Gly Gln Val Glu Ile Phe Pro Glu Ser Asp 560 565 570

GAT TTA GTG CAA TCC ATT AAT TTT GGT GAT AAT AGT GTT TAC CTG CCA 3750  
 Asp Leu Val Gln Ser Ile Asn Phe Gly Asp Asn Ser Val Tyr Leu Pro 575 580 585

ATA TTG AAT GAC TCT CAT GTA AAA AAC ATT ATT GAT TGT AAT GGA AAT 3798  
 Ile Leu Asn Asp Ser His Val Lys Asn Ile Ile Asp Cys Asn Gly Asn 590 595 600 605

AAC GAA TTA CGG TTG CAT AAC ATT GTC AAT TTT CTC TAT ACG GTA GGG 3846  
 Asn Glu Leu Arg Leu His Asn Ile Val Asn Phe Leu Tyr Thr Val Gly 610 615 620

CAA AGA TGG AAA AAT GAA GAA TAC TCA AGA CGC AGG ACA TAC ATT CGT 3894  
 Gln Arg Trp Lys Asn Glu Glu Tyr Ser Arg Arg Arg Thr Tyr Ile Arg 625 630 635

GAC TTA AAA AAA TAT ATG GGA TAT TCA GAA GAA ATG GCT AAG CTA GAG 3942  
 Asp Leu Lys Lys Tyr Met Gly Tyr Ser Glu Glu Met Ala Lys Leu Glu 640 645 650

GCC AAT TGG ATA TCT ATG ATT TTA TGT TCT AAA GGC GGC CTT TAT GAT 3990  
 Ala Asn Trp Ile Ser Met Ile Leu Cys Ser Lys Gly Gly Leu Tyr Asp 655 660 665

GTT GTA GAA AAT GAA CTT GGT TCT CGC CAT ATC ATG GAT GAA TGG CTA 4038  
 Val Val Glu Asn Glu Leu Gly Ser Arg His Ile Met Asp Glu Trp Leu 670 675 680 685

CCT CAG GAT GAA AGT TAT GTT CGG GCT TTT CCG AAA GGT AAA TCT GTA 4086  
 Pro Gln Asp Glu Ser Tyr Val Arg Ala Phe Pro Lys Gly Lys Ser Val 690 695 700

CAT CTG TTG GCA GGT AAT GTT CCA TTA TCT GGG ATC ATG TCT ATA TTA 4134  
 His Leu Leu Ala Gly Asn Val Pro Leu Ser Gly Ile Met Ser Ile Leu 705 710 715

CGC GCA ATT TTA ACT AAG AAT CAG TGT ATT ATA AAA ACA TCG TCA ACC 4182  
 Arg Ala Ile Leu Thr Lys Asn Gln Cys Ile Ile Lys Thr Ser Ser Thr 720 725 730

GAT CCT TTT ACC GCT AAT GCA TTA GCG TTA AGT TTT ATT GAT GTA GAC 4230  
 Asp Pro Phe Thr Ala Asn Ala Leu Ala Leu Ser Phe Ile Asp Val Asp 735 740 745

CCT AAT CAT CCG ATA ACG CGC TCT TTA TCT GTT ATA TAT TGG CCC CAC 4278  
 Pro Asn His Pro Ile Thr Arg Ser Leu Ser Val Ile Tyr Trp Pro His 750 755 760 765

CAA GGT GAT ACA TCA CTC GCA AAA GAA ATT ATG CGA CAT GCG GAT GTT 4326  
 Gln Gly Asp Thr Ser Leu Ala Lys Glu Ile Met Arg His Ala Asp Val 770 775 780

ATT GTC GCT TGG GGA GGG CCA GAT GCG ATT AAT TGG GCG GTA GAG CAT 4374  
 Ile Val Ala Trp Gly Gly Pro Asp Ala Ile Asn Trp Ala Val Glu His 785 790 795

GCG CCA TCT TAT GCT GAT GTG ATT AAA TTT GGT TCT AAA AAG AGT CTT 4422  
 Ala Pro Ser Tyr Ala Asp Val Ile Lys Phe Gly Ser Lys Lys Ser Leu 800 805 810

004240 29662960

TGC ATT ATC GAT AAT CCT GTT GAT TTG ACG TCC GCA GCG ACA GGT GCG 4470  
 Cys Ile Ile Asp Asn Pro Val Asp Leu Thr Ser Ala Ala Thr Gly Ala  
 815 820 825

GCT CAT GAT GTT TGT TTT TAC GAT CAG CGA GCT TGT TTT TCT GCC CAA 4518  
 Ala His Asp Val Cys Phe Tyr Asp Gln Arg Ala Cys Phe Ser Ala Gln  
 830 835 840 845

AAC ATA TAT TAC ATG GGA AAT CAT TAT GAG GAA TTT AAG TTA GCG TTG 4566  
 Asn Ile Tyr Tyr Met Gly Asn His Tyr Glu Glu Phe Lys Leu Ala Leu  
 850 855 860

ATA GAA AAA CTT AAT CTA TAT GCG CAT ATA TTA CCG AAT GCC AAA AAA 4614  
 Ile Glu Lys Leu Asn Leu Tyr Ala His Ile Leu Pro Asn Ala Lys Lys  
 865 870 875

GAT TTT GAT GAA AAG GCG GCC TAT TCT TTA GTT CAA AAA GAA AGC TTG 4662  
 Asp Phe Asp Glu Lys Ala Ala Tyr Ser Leu Val Gln Lys Glu Ser Leu  
 880 885 890

TTT GCT GGA TTA AAA GTA GAG GTG GAT ATT CAT CAA CGT TGG ATG ATT 4710  
 Phe Ala Gly Leu Lys Val Glu Val Asp Ile His Gln Arg Trp Met Ile  
 895 900 905

ATT GAG TCA AAT GCA GGT GTG GAA TTT AAT CAA CCA CTT GGC AGA TGT 4758  
 Ile Glu Ser Asn Ala Gly Val Glu Phe Asn Gln Pro Leu Gly Arg Cys  
 910 915 920 925

GTG TAC CTT CAT CAC GTC GAT AAT ATT GAG CAA ATA TTG CCT TAT GTT 4806  
 Val Tyr Leu His His Val Asp Asn Ile Glu Gln Ile Leu Pro Tyr Val  
 930 935 940

CAA AAA AAT AAG ACG CAA ACC ATA TCT ATT TTT CCT TGG GAG TCA TCA 4854  
 Gln Lys Asn Lys Thr Gln Thr Ile Ser Ile Phe Pro Trp Glu Ser Ser  
 945 950 955

TTT AAA TAT CGA GAT GCG TTA GCA TTA AAA GGT GCG GAA AGG ATT GTA 4902  
 Phe Lys Tyr Arg Asp Ala Leu Ala Leu Lys Gly Ala Glu Arg Ile Val  
 960 965 970

GAA GCA GGA ATG AAT AAC ATA TTT CGA GTT GGT GGA TCT CAT GAC GGA 4950  
 Glu Ala Gly Met Asn Asn Ile Phe Arg Val Gly Gly Ser His Asp Gly  
 975 980 985

ATG AGA CCG TTG CAA CGA TTA GTG ACA TAT ATT TCT CAT GAA AGG CCA 4998  
 Met Arg Pro Leu Gln Arg Leu Val Thr Tyr Ile Ser His Glu Arg Pro  
 990 995 1000 1005

TCT AAC TAT ACG GCT AAG GAT GTT GCG GTT GAA ATA GAA CAG ACT CGA 5046  
 Ser Asn Tyr Thr Ala Lys Asp Val Ala Val Glu Ile Glu Gln Thr Arg  
 1010 1015 1020

TTC CTG GAA GAA GAT AAG TTC CTT GTA TTT GTC CCA TAATAGGTAA 5092  
 Phe Leu Glu Glu Asp Lys Phe Leu Val Phe Val Pro  
 1025 1030

AAGT ATG GAA AAT GAA TCA AAA TAT AAA ACC ATC GAC CAC GTT ATT TGT 5141  
 Met Glu Asn Glu Ser Lys Tyr Lys Thr Ile Asp His Val Ile Cys  
 1 5 10 15

GTT GAA GGA AAT AAA AAA ATT CAT GTT TGG GAA ACG CTG CCA GAA GAA 5189  
 Val Glu Gly Asn Lys Lys Ile His Val Trp Glu Thr Leu Pro Glu Glu  
 20 25 30

AAC AGC CCA AAG AGA AAG AAT GCC ATT ATT ATT GCG TCT GGT TTT GCC 5237  
 Asn Ser Pro Lys Arg Lys Asn Ala Ile Ile Ile Ala Ser Gly Phe Ala  
 35 40 45

004240 29662560

CGC AGG ATG GAT CAT TTT GCT GGT CTG GCG GAA TAT TTA TCG CGG AAT 5285  
 Arg Arg Met Asp His Phe Ala Gly Leu Ala Glu Tyr Leu Ser Arg Asn  
 50 55 60

GGA TTT CAT GTG ATC CGC TAT GAT TCG CTT CAC CAC GTT GGA TTG AGT 5333  
 Gly Phe His Val Ile Arg Tyr Asp Ser Leu His His Val Gly Leu Ser  
 65 70 75

TCA GGG ACA ATT GAT GAA TTT ACA ATG TCT ATA GGA AAG CAG AGC TTG 5381  
 Ser Gly Thr Ile Asp Glu Phe Thr Met Ser Ile Gly Lys Gln Ser Leu  
 80 85 90 95

TTA GCA GTG GTT GAT TGG TTA ACT ACA CGA AAA ATA AAT AAC TTC GGT 5429  
 Leu Ala Val Val Asp Trp Leu Thr Thr Arg Lys Ile Asn Asn Phe Gly  
 100 105 110

ATG TTG GCT TCA AGC TTA TCT GCG CGG ATA GCT TAT GCA AGC CTA TCT 5477  
 Met Leu Ala Ser Ser Leu Ser Ala Arg Ile Ala Tyr Ala Ser Leu Ser  
 115 120 125

GAA ATC AAT GCT TCG TTT TTA ATC ACC GCA GTC GGT GTT GTT AAC TTA 5525  
 Glu Ile Asn Ala Ser Phe Leu Ile Thr Ala Val Gly Val Val Asn Leu  
 130 135 140

AGA TAT TCT CTT GAA AGA GCT TTA GGG TTT GAT TAT CTC AGT CTA CCC 5573  
 Arg Tyr Ser Leu Glu Arg Ala Leu Gly Phe Asp Tyr Leu Ser Leu Pro  
 145 150 155

ATT AAT GAA TTG CCG GAT AAT CTA GAT TTT GAA GGC CAT AAA TTG GGT 5621  
 Ile Asn Glu Leu Pro Asp Asn Leu Asp Phe Glu Gly His Lys Leu Gly  
 160 165 170 175

GCT GAA GTC TTT GCG AGA GAT TGT CTT GAT TTT GGT TGG GAA GAT TTA 5669  
 Ala Glu Val Phe Ala Arg Asp Cys Leu Asp Phe Gly Trp Glu Asp Leu  
 180 185 190

GCT TCT ACA ATT AAT AAC ATG ATG TAT CTT GAT ATA CCG TTT ATT GCT 5717  
 Ala Ser Thr Ile Asn Asn Met Met Tyr Leu Asp Ile Pro Phe Ile Ala  
 195 200 205

TTT ACT GCA AAT AAC GAT AAT TGG GTC AAG CAA GAT GAA GTT ATC ACA 5765  
 Phe Thr Ala Asn Asn Asp Asn Trp Val Lys Gln Asp Glu Val Ile Thr  
 210 215 220

TTG TTA TCA AAT ATT CGT AGT AAT CGA TGC AAG ATA TAT TCT TTG TTA 5813  
 Leu Leu Ser Asn Ile Arg Ser Asn Arg Cys Lys Ile Tyr Ser Leu Leu  
 225 230 235

GGA AGT TCG CAT GAC TTG AGT GAA AAT TTA GTG GTC CTG CGC AAT TTT 5861  
 Gly Ser Ser His Asp Leu Ser Glu Asn Leu Val Val Leu Arg Asn Phe  
 240 245 250 255

TAT CAA TCG GTT ACG AAA GCC GCT ATC GCG ATG GAT AAT GAT CAT CTG 5909  
 Tyr Gln Ser Val Thr Lys Ala Ala Ile Ala Met Asp Asn Asp His Leu  
 260 265 270

GAT ATT GAT GTT GAT ATT ACT GAA CCG TCA TTT GAA CAT TTA ACT ATT 5957  
 Asp Ile Asp Val Asp Ile Thr Glu Pro Ser Phe Glu His Leu Thr Ile  
 275 280 285

GCG ACA GTC AAT GAA CGC CGA ATG AGA ATT GAG ATT GAA AAT CAA GCA 6005  
 Ala Thr Val Asn Glu Arg Met Arg Ile Glu Ile Glu Asn Gln Ala  
 290 295 300

ATT TCT CTG TCT TAAATCTAT TGAGATATTC TATCACTCAA ATAGCAATAT 6057  
 Ile Ser Leu Ser  
 305

004240" 29662560



AAGGACTCTC	T	ATG	AAA	TTT	GGA	AAC	TTT	TTG	CTT	ACA	TAC	CAA	CCT	CCC	6107	
		Met	Lys	Phe	Gly	Asn	Phe	Leu	Leu	Thr	Tyr	Gln	Pro	Pro		
		1				5					10					
CAA	TTT	TCT	CAA	ACA	GAG	GTA	ATG	AAA	CGT	TTG	GTT	AAA	TTA	GCT	CGC	6155
Gln	Phe	Ser	Gln	Thr	Glu	Val	Met	Lys	Arg	Leu	Val	Lys	Leu	Gly	Arg	
	15					20					25					
ATC	TCT	GAG	GAG	TGT	GGT	TTT	GAT	ACC	GTA	TGG	TTA	CTG	GAG	CAT	CAT	6203
Ile	Ser	Glu	Glu	Cys	Gly	Phe	Asp	Thr	Val	Trp	Leu	Leu	Glu	His	His	
	30				35					40					45	
TTC	ACG	GAG	TTT	GGT	TTG	CTT	GGT	AAC	CCT	TAT	GTC	GCT	GCT	GCA	TAT	6251
Phe	Thr	Glu	Phe	Gly	Leu	Leu	Gly	Asn	Pro	Tyr	Val	Ala	Ala	Ala	Tyr	
				50					55					60		
TTA	CTT	GGC	GCG	ACT	AAA	AAA	TTG	AAT	GTA	GGA	ACT	GCC	GCT	ATT	GTT	6299
Leu	Leu	Gly	Ala	Thr	Lys	Lys	Leu	Asn	Val	Gly	Thr	Ala	Ala	Ile	Val	
			65					70					75			
CTT	CCC	ACA	GCC	CAT	CCA	GTA	CGC	CAA	CTT	GAA	GAT	GTG	AAT	TTA	TTG	6347
Leu	Pro	Thr	Ala	His	Pro	Val	Arg	Gln	Leu	Glu	Asp	Val	Asn	Leu	Leu	
		80					85					90				
GAT	CAA	ATG	TCA	AAA	GGA	CGA	TTT	CGG	TTT	GGT	ATT	TGC	CGA	GGG	CTT	6395
Asp	Gln	Met	Ser	Lys	Gly	Arg	Phe	Arg	Phe	Gly	Ile	Cys	Arg	Gly	Leu	
	95					100					105					
TAC	AAC	AAG	GAC	TTT	CGC	GTA	TTT	GGC	ACA	GAT	ATG	AAT	AAC	AGT	CGC	6443
Tyr	Asn	Lys	Asp	Phe	Arg	Val	Phe	Gly	Thr	Asp	Met	Asn	Asn	Ser	Arg	
	110				115					120					125	
GCC	TTA	GCG	GAA	TGC	TGG	TAC	GGG	CTG	ATA	AAG	AAT	GGC	ATG	ACA	GAG	6491
Ala	Leu	Ala	Glu	Cys	Trp	Tyr	Gly	Leu	Ile	Lys	Asn	Gly	Met	Thr	Glu	
				130					135					140		
GGA	TAT	ATG	GAA	GCT	GAT	AAT	GAA	CAT	ATC	AAG	TTC	CAT	AAG	GTA	AAA	6539
Gly	Tyr	Met	Glu	Ala	Asp	Asn	Glu	His	Ile	Lys	Phe	His	Lys	Val	Lys	
			145					150					155			
GTA	AAC	CCC	GCG	GCG	TAT	AGC	AGA	GGT	GGC	GCA	CCG	GTT	TAT	GTG	GTG	6587
Val	Asn	Pro	Ala	Ala	Tyr	Ser	Arg	Gly	Gly	Ala	Pro	Val	Tyr	Val	Val	
		160					165					170				
GCT	GAA	TCA	GCT	TGC	ACG	ACT	GAG	TGG	GCT	GCT	CAA	TTT	GGC	CTA	CCG	6635
Ala	Glu	Ser	Ala	Ser	Thr	Thr	Glu	Trp	Ala	Ala	Gln	Phe	Gly	Leu	Pro	
	175					180					185					
ATG	ATA	TTA	AGT	TGG	ATT	ATA	AAT	ACT	AAC	GAA	AAG	AAA	GCA	CAA	CTT	6683
Met	Ile	Leu	Ser	Trp	Ile	Ile	Asn	Thr	Asn	Glu	Lys	Lys	Ala	Gln	Leu	
	190				195				200						205	
GAG	CTT	TAT	AAT	GAA	GTG	GCT										

[illegible]

TAT GAT TTC AAT AAA GGG CAG TGG CGT GAC TTT GTA TTA AAA GGA CAT 6923  
 Tyr Asp Phe Asn Lys Gly Gln Trp Arg Asp Phe Val Leu Lys Gly His 285  
 270 275 280

AAA GAT ACT AAT CGC CGT ATT GAT TAC AGT TAC GAA ATC AAT CCC GTG 6971  
 Lys Asp Thr Asn Arg Arg Ile Asp Tyr Ser Tyr Glu Ile Asn Pro Val 300  
 290 295

GGA ACG CCG CAG GAA TGT ATT GAC ATA ATT CAA AAA GAC ATT GAT GCT 7019  
 Gly Thr Pro Gln Glu Cys Ile Asp Ile Ile Gln Lys Asp Ile Asp Ala 315  
 305 310

ACA GGA ATA TCA AAT ATT TGT TGT GGA TTT GAA GCT AAT GGA ACA GTA 7067  
 Thr Gly Ile Ser Asn Ile Cys Cys Gly Phe Glu Ala Asn Gly Thr Val 330  
 320 325

GAC GAA ATT ATT GCT TCC ATG AAG CTC TTC CAG TCT GAT GTC ATG CCA 7115  
 Asp Glu Ile Ile Ala Ser Met Lys Leu Phe Gln Ser Asp Val Met Pro 345  
 335 340

TTT CTT AAA GAA AAA CAA CGT TCG CTA TTA TAT TAGCTAAGGA GAAAGAA 7165  
 Phe Leu Lys Glu Lys Gln Arg Ser Leu Leu Tyr 360  
 350 355

ATG AAA TTT GGA TTG TTC TTC CTT AAC TTC ATC AAT TCA ACA ACT GTT 7213  
 Met Lys Phe Gly Leu Phe Phe Leu Asn Phe Ile Asn Ser Thr Thr Val 15  
 1 5 10

CAA GAA CAA AGT ATA GTT CGC ATG CAG GAA ATA ACG GAG TAT GTT GAT 7261  
 Gln Glu Gln Ser Ile Val Arg Met Gln Glu Ile Thr Glu Tyr Val Asp 30  
 20 25

AAG TTG AAT TTT GAA CAG ATT TTA GTG TAT GAA AAT CAT TTT TCA GAT 7309  
 Lys Leu Asn Phe Glu Gln Ile Leu Val Tyr Glu Asn His Phe Ser Asp 45  
 35 40

AAT GGT GTT GTC GGC GCT CCT CTG ACT GTT TCT GGT TTT CTG CTC GGT 7357  
 Asn Gly Val Val Gly Ala Pro Leu Thr Val Ser Gly Phe Leu Leu Gly 60  
 50 55

TTA ACA GAG AAA ATT AAA ATT GGT TCA TTA AAT CAC ATC ATT ACA ACT 7405  
 Leu Thr Glu Lys Ile Lys Ile Gly Ser Leu Asn His Ile Ile Thr Thr 80  
 65 70 75

CAT CAT CCT GTC GCC ATA GCG GAG GAA GCT TGC TTA TTG GAT CAG TTA 7453  
 His His Pro Val Ala Ile Ala Glu Glu Ala Cys Leu Leu Asp Gln Leu 95  
 85 90

AGT GAA GGG AGA TTT ATT TTA GGG TTT AGT GAT TGC GAA AAA AAA GAT 7501  
 Ser Glu Gly Arg Phe Ile Leu Gly Phe Ser Asp Cys Glu Lys Lys Asp 110  
 100 105

GAA ATG CAT TTT TTT AAT CGC CCG GTT GAA TAT CAA CAG CAA CTA TTT 7549  
 Glu Met His Phe Phe Asn Arg Pro Val Glu Tyr Gln Gln Gln Leu Phe 125  
 115 120

GAA GAG TGT TAT GAA ATC ATT AAC GAT GCT TTA ACA ACA GGC TAT TGT 7597  
 Glu Glu Cys Tyr Glu Ile Ile Asn Asp Ala Leu Thr Thr Gly Tyr Cys 140  
 130 135

AAT CCA GAT AAC GAT TTT TAT AGC TTC CCT AAA ATA TCT GTA AAT CCC 7645  
 Asn Pro Asp Asn Asp Phe Tyr Ser Phe Pro Lys Ile Ser Val Asn Pro 160  
 145 150 155

CAT GCT TAT ACG CCA GGC GGA CCT CGG AAA TAT GTA ACA GCA ACC AGT 7693  
 His Ala Tyr Thr Pro Gly Gly Pro Arg Lys Tyr Val Thr Ala Thr Ser 175  
 165 170

004240 0962560

CAT CAT ATT GTT GAG TGG GCG GCC AAA AAA GGT ATT CCT CTC ATC TTT 7741  
His His Ile Val Glu Trp Ala Ala Lys Lys Gly Ile Pro Leu Ile Phe  
180 185 190

AAG TGG GAT GAT TCT AAT GAT GTT AGA TAT GAA TAT GCT GAA AGA TAT 7789  
Lys Trp Asp Asp Ser Asn Asp Val Arg Tyr Glu Tyr Ala Glu Arg Tyr  
195 200 205

AAA GCC GTT GCG GAT AAA TAT GAC GTT GAC CTA TCA GAG ATA GAC CAT 7837  
Lys Ala Val Ala Asp Lys Tyr Asp Val Asp Leu Ser Glu Ile Asp His  
210 215 220

CAG TTA ATG ATA TTA GTT AAC TAT AAC GAA GAT AGT AAT AAA GCT AAA 7885  
Gln Leu Met Ile Leu Val Asn Tyr Asn Glu Asp Ser Asn Lys Ala Lys  
225 230 235 240

CAA GAG ACG CGT GCA TTT ATT AGT GAT TAT GTT CTT GAA ATG CAC CCT 7933  
Gln Glu Thr Arg Ala Phe Ile Ser Asp Tyr Val Leu Glu Met His Pro  
245 250 255

AAT GAA AAT TTC GAA AAT AAA CTT GAA GAA ATA ATT GCA GAA AAC GCT 7981  
Asn Glu Asn Phe Glu Asn Lys Leu Glu Glu Ile Ile Ala Glu Asn Ala  
260 265 270

GTC GGA AAT TAT ACG GAG TGT ATA ACT GCG GCT AAG TTG GCA ATT GAA 8029  
Val Gly Asn Tyr Thr Glu Cys Ile Thr Ala Ala Lys Leu Ala Ile Glu  
275 280 285

AAG TGT GGT GCG AAA AGT GTA TTG CTG TCC TTT GAA CCA ATG AAT GAT 8077  
Lys Cys Gly Ala Lys Ser Val Leu Leu Ser Phe Glu Pro Met Asn Asp  
290 295 300

TTG ATG AGC CAA AAA AAT GTA ATC AAT ATT GTT GAT GAT AAT ATT AAG 8125  
Leu Met Ser Gln Lys Asn Val Ile Asn Ile Val Asp Asp Asn Ile Lys  
305 310 315 320

AAG TAC CAC ATG GAA TAT ACC TAATAGATTT CGAGTTGCAG CGAGGCCGCA 8176  
Lys Tyr His Met Glu Tyr Thr  
325

AGTGAACGAA TCCOCAGGAG CATAGATAAC TATGTGACTG GGGTGAGTGA AAGCAGCCAA 8236

CAAAGCAGCA GCTTGAAAG ATG AAG GGT ATA AAA GAG TAT GAC AGC AGT GCT 8288  
Met Lys Gly Ile Lys Glu Tyr Asp Ser Ser Ala  
1 5 10

GCC ATA CTT TCT AAT ATT ATC TTG AGG AGT AAA ACA GGT ATG ACT TCA 8336  
Ala Ile Leu Ser Asn Ile Ile Leu Arg Ser Lys Thr Gly Met Thr Ser  
15 20 25

TAT GTT GAT AAA CAA GAA ATT ACA GCA AGC TCA GAA ATT GAT GAT TTG 8384  
Tyr Val Asp Lys Gln Glu Ile Thr Ala Ser Ser Glu Ile Asp Asp Leu  
30 35 40

ATT TTT TCG AGC GAT CCA TTA GTG TGG TCT TAC GAC GAG CAG GAA AAA 8432  
Ile Phe Ser Ser Asp Pro Leu Val Trp Ser Tyr Asp Glu Gln Glu Lys  
45 50 55

ATC AGA AAG AAA CTT GTG CTT GAT GCA TTT CGT AAT CAT TAT AAA CAT 8480  
Ile Arg Lys Lys Leu Val Leu Asp Ala Phe Arg Asn His Tyr Lys His  
60 65 70 75

TGT CGA GAA TAT CGT CAC TAC TGT CAG GCA CAC AAA GTA GAT GAC AAT 8528  
Cys Arg Glu Tyr Arg His Tyr Cys Gln Ala His Lys Val Asp Asp Asn  
80 85 90

00440" 29662560

ATT ACG GAA ATT GAT GAC ATA CCT GTA TTC CCA ACA TCG GTT TTT AAG 8576  
 Ile Thr Glu Ile Asp Asp Ile Pro Val Phe Pro Thr Ser Val Phe Lys  
 95 100 105

TTT ACT CGC TTA TTA ACT TCT CAG GAA AAC GAG ATT GAA AGT TGG TTT 8624  
 Phe Thr Arg Leu Leu Thr Ser Gln Glu Asn Glu Ile Glu Ser Trp Phe  
 110 115 120

ACC AGT AGC GGC ACG AAT GGT TTA AAA AGT CAG GTG GCG CGT GAC AGA 8672  
 Thr Ser Ser Gly Thr Asn Gly Leu Lys Ser Gln Val Ala Arg Asp Arg  
 125 130 135

TTA AGT ATT GAG AGA CTC TTA GGC TCT GTG AGT TAT GGC ATG AAA TAT 8720  
 Leu Ser Ile Glu Arg Leu Leu Gly Ser Val Ser Tyr Gly Met Lys Tyr  
 140 145 150 155

GTT GGT AGT TGG TTT GAT CAT CAA ATA GAA TTA GTC AAT TTG GGA CCA 8768  
 Val Gly Ser Trp Phe Asp His Gln Ile Glu Leu Val Asn Leu Gly Pro  
 160 165 170

GAT AGA TTT AAT GCT CAT AAT ATT TGG TTT AAA TAT GTT ATG AGT TTG 8816  
 Asp Arg Phe Asn Ala His Asn Ile Trp Phe Lys Tyr Val Met Ser Leu  
 175 180 185

GTG GAA TTG TTA TAT CCT ACG ACA TTT ACC GTA ACA GAA GAA CGA ATA 8864  
 Val Glu Leu Leu Tyr Pro Thr Phe Thr Val Thr Glu Glu Arg Ile  
 190 195 200

GAT TTT GTT AAA ACA TTG AAT AGT CTT GAA CGA ATA AAA AAT CAA GGG 8912  
 Asp Phe Val Lys Thr Leu Asn Ser Leu Glu Arg Ile Lys Asn Gln Gly  
 205 210 215

AAA GAT CTT TGT CTT ATT GGT TCG CCA TAC TTT ATT TAT TTA CTC TGC 8960  
 Lys Asp Leu Cys Leu Ile Gly Ser Pro Tyr Phe Ile Tyr Leu Leu Cys  
 220 225 230 235

CAT TAT ATG AAA GAT AAA AAA ATC TCA TTT TCT GGA GAT AAA AGC CTT 9008  
 His Tyr Met Lys Asp Lys Lys Ile Ser Phe Ser Gly Asp Lys Ser Leu  
 240 245 250

TAT ATC ATA ACC GCA GGC GGC TGG AAA AGT TAC GAA AAA GAA TCT CTG 9056  
 Tyr Ile Ile Thr Gly Gly Gly Trp Lys Ser Tyr Glu Lys Glu Ser Leu  
 255 260 265

AAA CGT GAT GAT TTC AAT CAT CTT TTA TTT GAT ACT TTC AAT CTC AGT 9104  
 Lys Arg Asp Asp Phe Asn His Leu Leu Phe Asp Thr Phe Asn Leu Ser  
 270 275 280

GAT ATT AGT CAG ATC CGA GAT ATA TTT AAT CAA GTT GAA CTC AAC ACT 9152  
 Asp Ile Ser Gln Ile Arg Asp Ile Phe Asn Gln Val Glu Leu Asn Thr  
 285 290 295

TGT TTC TTT GAG GAT GAA ATG CAG CGT AAA CAT GTT CCG CCG TGG GTA 9200  
 Cys Phe Phe Glu Asp Glu Met Gln Arg Lys His Val Pro Pro Trp Val  
 300 305 310 315

TAT GCG CGA GCG CTT GAT CCT GAA ACG TTG AAA CCT GTA CCT GAT GGA 9248  
 Tyr Ala Arg Ala Leu Asp Pro Glu Thr Leu Lys Pro Val Pro Asp Gly  
 320 325 330

ACG CCG GGG TTG ATG AGT TAT ATG GAT GCG TCA GCA ACC AGT TAT CCA 9296  
 Thr Pro Gly Leu Met Ser Tyr Met Asp Ala Ser Ala Thr Ser Tyr Pro  
 335 340 345

GCA TTT ATT GTT ACC GAT GAT GTC GGG ATA ATT AGC AGA GAA TAT GGT 9344  
 Ala Phe Ile Val Thr Asp Asp Val Gly Ile Ile Ser Arg Glu Tyr Gly  
 350 355 360

004240 29662560

AAG TAT CCC GGC GTG CTC GTT GAA ATT TTA CGT CGC GTC AAT ACG AGG 9392  
 Lys Tyr Pro Gly Val Leu Val Glu Ile Leu Arg Arg Val Asn Thr Arg  
 365 370 375  
 ACG CAG AAA GGG TGT GCT TTA AGC TTA ACC GAA GCG TTT GAT AGT 9437  
 Thr Gln Lys Gly Cys Ala Leu Ser Leu Thr Glu Ala Phe Asp Ser  
 380 385 390  
 TGATATCCTT TGCCTAATTG TAAGTGAAT GCTTGCCTTA TATAAATCTG AATGACATCT 9497  
 ACACTTTACA AAATTCTCCA AAACATCCAC ATTTGGGTAC TTGATAGAGG TTTATGGGGT 9557  
 TGGCTTAACA TTGTTCTCAT TGTTATTATT GGCTCAAAGC AAAAGGAGAT AACATGAAAA 9617  
 AATTGGCAGT TATGCTTGCA TTGGGAATGA TTAGCTTTGG TCAATGGCA GTTGATGGGT 9677  
 ATAAAGATGC AAAGTTTGGC ATGACAGAAG AAGAGTTTCT TTCGAAGAGG TTATGTGATT 9737  
 TTGAAAAATT TGAGGGAGAT TCTCGAATAG AAGAAGTATC ACTTTATTCA TGTTCGACT 9797  
 TTTCGTTTGC TAACAAAAAG CGTGAAGCAA TGGCATTTTT TTAAATGGG AAATTTAAAA 9857  
 GATTAGAGAT TAATATTGGC AGACTTGTA AGCCAGTAAG CAAATCGTTA ACGAAAAAGT 9917  
 ACGGAGATGG ATCATCGTAT CCATCAAAAG AAGAATTTGA GAACGCGCTA AAATACAATG 9977  
 GAACTATGTC TATAGGTTAT GATAATAATA CGGTATTAGT TGATATACAT ATAATATGTG 10037  
 GCAAAGAAGG CATAGAAACC AGTCAACTGA TTTATACGAG TCCAGATGTT TATACGCTCC 10097  
 CAGATTTCGG AGAAAAATC CAGGAATTAA AGGGATTAAA GGAATTCGAG CTCGGTACCC 10157  
 GGGGATCCCT CGAGGTCGAC CTGCAGGCAG CECTTGGCGT CACCCGCAGT TCGGTGGTTA 10217  
 ATA 10220

## (2) INFORMATION FOR SEQ ID NO: 4:

- (i) SEQUENCE CHARACTERISTICS:
- (A) LENGTH: 483 amino acids
  - (B) TYPE: amino acid
  - (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 4:

Met Ala Asn Met Thr Lys Lys Ile Ser Phe Ile Ile Asn Gly Gln Val  
 1 5 10 15  
 Glu Ile Phe Pro Glu Ser Asp Asp Leu Val Gln Ser Ile Asn Phe Gly  
 20 25 30  
 Asp Asn Ser Val Tyr Leu Pro Ile Leu Asn Asp Ser His Val Lys Asn  
 35 40 45  
 Ile Ile Asp Cys Asn Gly Asn Asn Glu Leu Arg Leu His Asn Ile Val  
 50 55 60  
 Asn Phe Leu Tyr Thr Val Gly Gln Arg Trp Lys Asn Glu Glu Tyr Ser  
 65 70 75 80  
 Arg Arg Arg Thr Tyr Ile Arg Asp Leu Lys Lys Tyr Met Gly Tyr Ser  
 85 90 95  
 Glu Glu Met Ala Lys Leu Glu Ala Asn Trp Ile Ser Met Ile Leu Cys  
 100 105 110

004240" 29662660

Ser Lys Gly Gly Leu Tyr Asp Val Val Glu Asn Glu Leu Gly Ser Arg  
 115 120 125  
 His Ile Met Asp Glu Trp Leu Pro Gln Asp Glu Ser Tyr Val Arg Ala  
 130 135 140  
 Phe Pro Lys Gly Lys Ser Val His Leu Leu Ala Gly Asn Val Pro Leu  
 145 150 155 160  
 Ser Gly Ile Met Ser Ile Leu Arg Ala Ile Leu Thr Lys Asn Gln Cys  
 165 170 175  
 Ile Ile Lys Thr Ser Ser Thr Asp Pro Phe Thr Ala Asn Ala Leu Ala  
 180 185 190  
 Leu Ser Phe Ile Asp Val Asp Pro Asn His Pro Ile Thr Arg Ser Leu  
 195 200 205  
 Ser Val Ile Tyr Trp Pro His Gln Gly Asp Thr Ser Leu Ala Lys Glu  
 210 215 220  
 Ile Met Arg His Ala Asp Val Ile Val Ala Trp Gly Gly Pro Asp Ala  
 225 230 235 240  
 Ile Asn Trp Ala Val Glu His Ala Pro Ser Tyr Ala Asp Val Ile Lys  
 245 250 255  
 Phe Gly Ser Lys Lys Ser Leu Cys Ile Ile Asp Asn Pro Val Asp Leu  
 260 265 270  
 Thr Ser Ala Ala Thr Gly Ala Ala His Asp Val Cys Phe Tyr Asp Gln  
 275 280 285  
 Arg Ala Cys Phe Ser Ala Gln Asn Ile Tyr Tyr Met Gly Asn His Tyr  
 290 295 300  
 Glu Glu Phe Lys Leu Ala Leu Ile Glu Lys Leu Asn Leu Tyr Ala His  
 305 310 315 320  
 Ile Leu Pro Asn Ala Lys Lys Asp Phe Asp Glu Lys Ala Ala Tyr Ser  
 325 330 335  
 Leu Val Gln Lys Glu Ser Leu Phe Ala Gly Leu Lys Val Glu Val Asp  
 340 345 350  
 Ile His Gln Arg Trp Met Ile Ile Glu Ser Asn Ala Gly Val Glu Phe  
 355 360 365  
 Asn Gln Pro Leu Gly Arg Cys Val Tyr Leu His His Val Asp Asn Ile  
 370 375 380  
 Glu Gln Ile Leu Pro Tyr Val Gln Lys Asn Lys Thr Gln Thr Ile Ser  
 385 390 395 400  
 Ile Phe Pro Trp Glu Ser Ser Phe Lys Tyr Arg Asp Ala Leu Ala Leu  
 405 410 415  
 Lys Gly Ala Glu Arg Ile Val Glu Ala Gly Met Asn Asn Ile Phe Arg  
 420 425 430  
 Val Gly Gly Ser His Asp Gly Met Arg Pro Leu Gln Arg Leu Val Thr  
 435 440 445  
 Tyr Ile Ser His Glu Arg Pro Ser Asn Tyr Thr Ala Lys Asp Val Ala  
 450 455 460  
 Val Glu Ile Glu Gln Thr Arg Phe Leu Glu Glu Asp Lys Phe Leu Val  
 465 470 475 480

Phe Val Pro

## (2) INFORMATION FOR SEQ ID NO: 5:

## (i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 307 amino acids
- (B) TYPE: amino acid
- (D) TOPOLOGY: linear

## (ii) MOLECULE TYPE: protein

## (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 5:

Met Glu Asn Glu Ser Lys Tyr Lys Thr Ile Asp His Val Ile Cys Val  
 1 5 10 15  
 Glu Gly Asn Lys Lys Ile His Val Trp Glu Thr Leu Pro Glu Glu Asn  
 20 25 30  
 Ser Pro Lys Arg Lys Asn Ala Ile Ile Ile Ala Ser Gly Phe Ala Arg  
 35 40 45  
 Arg Met Asp His Phe Ala Gly Leu Ala Glu Tyr Leu Ser Arg Asn Gly  
 50 55 60  
 Phe His Val Ile Arg Tyr Asp Ser Leu His His Val Gly Leu Ser Ser  
 65 70 75 80  
 Gly Thr Ile Asp Glu Phe Thr Met Ser Ile Gly Lys Gln Ser Leu Leu  
 85 90 95  
 Ala Val Val Asp Trp Leu Thr Thr Arg Lys Ile Asn Asn Phe Gly Met  
 100 105 110  
 Leu Ala Ser Ser Leu Ser Ala Arg Ile Ala Tyr Ala Ser Leu Ser Glu  
 115 120 125  
 Ile Asn Ala Ser Phe Leu Ile Thr Ala Val Gly Val Val Asn Leu Arg  
 130 135 140  
 Tyr Ser Leu Glu Arg Ala Leu Gly Phe Asp Tyr Leu Ser Leu Pro Ile  
 145 150 155 160  
 Asn Glu Leu Pro Asp Asn Leu Asp Phe Glu Gly His Lys Leu Gly Ala  
 165 170 175  
 Glu Val Phe Ala Arg Asp Cys Leu Asp Phe Gly Trp Glu Asp Leu Ala  
 180 185 190  
 Ser Thr Ile Asn Asn Met Met Tyr Leu Asp Ile Pro Phe Ile Ala Phe  
 195 200 205  
 Thr Ala Asn Asn Asp Asn Trp Val Lys Gln Asp Glu Val Ile Thr Leu  
 210 215 220  
 Leu Ser Asn Ile Arg Ser Asn Arg Cys Lys Ile Tyr Ser Leu Leu Gly  
 225 230 235 240  
 Ser Ser His Asp Leu Ser Glu Asn Leu Val Val Leu Arg Asn Phe Tyr  
 245 250 255  
 Gln Ser Val Thr Lys Ala Ala Ile Ala Met Asp Asn Asp His Leu Asp  
 260 265 270  
 Ile Asp Val Asp Ile Thr Glu Pro Ser Phe Glu His Leu Thr Ile Ala  
 275 280 285

004240 29662560

Thr Val Asn Glu Arg Arg Met Arg Ile Glu Ile Glu Asn Gln Ala Ile  
 290 295 300

Ser Leu Ser  
 305

(2) INFORMATION FOR SEQ ID NO: 6:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 360 amino acids  
 (B) TYPE: amino acid  
 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 6:

Met Lys Phe Gly Asn Phe Leu Leu Thr Tyr Gln Pro Pro Gln Phe Ser  
 1 5 10 15  
 Gln Thr Glu Val Met Lys Arg Leu Val Lys Leu Gly Arg Ile Ser Glu  
 20 25 30  
 Glu Cys Gly Phe Asp Thr Val Trp Leu Leu Glu His His Phe Thr Glu  
 35 40 45  
 Phe Gly Leu Leu Gly Asn Pro Tyr Val Ala Ala Ala Tyr Leu Leu Gly  
 50 55 60  
 Ala Thr Lys Lys Leu Asn Val Gly Thr Ala Ala Ile Val Leu Pro Thr  
 65 70 75 80  
 Ala His Pro Val Arg Gln Leu Glu Asp Val Asn Leu Leu Asp Gln Met  
 85 90 95  
 Ser Lys Gly Arg Phe Arg Phe Gly Ile Cys Arg Gly Leu Tyr Asn Lys  
 100 105 110  
 Asp Phe Arg Val Phe Gly Thr Asp Met Asn Asn Ser Arg Ala Leu Ala  
 115 120 125  
 Glu Cys Trp Tyr Gly Leu Ile Lys Asn Gly Met Thr Glu Gly Tyr Met  
 130 135 140  
 Glu Ala Asp Asn Glu His Ile Lys Phe His Lys Val Lys Val Asn Pro  
 145 150 155 160  
 Ala Ala Tyr Ser Arg Gly Gly Ala Pro Val Tyr Val Val Ala Glu Ser  
 165 170 175  
 Ala Ser Thr Thr Glu Trp Ala Ala Gln Phe Gly Leu Pro Met Ile Leu  
 180 185 190  
 Ser Trp Ile Ile Asn Thr Asn Glu Lys Lys Ala Gln Leu Glu Leu Tyr  
 195 200 205  
 Asn Glu Val Ala Gln Glu Tyr Gly His Asp Ile His Asn Ile Asp His  
 210 215 220  
 Cys Leu Ser Tyr Ile Thr Ser Val Asp His Asp Ser Ile Lys Ala Lys  
 225 230 235 240  
 Glu Ile Cys Arg Lys Phe Leu Gly His Trp Tyr Asp Ser Tyr Val Asn  
 245 250 255  
 Ala Thr Thr Ile Phe Asp Asp Ser Asp Gln Thr Arg Gly Tyr Asp Phe  
 260 265 270

004240 49662560



Asn Lys Gly Gln Trp Arg Asp Phe Val Leu Lys Gly His Lys Asp Thr  
 275 280 285  
 Asn Arg Arg Ile Asp Tyr Ser Tyr Glu Ile Asn Pro Val Gly Thr Pro  
 290 295 300  
 Gln Glu Cys Ile Asp Ile Ile Gln Lys Asp Ile Asp Ala Thr Gly Ile  
 305 310 315 320  
 Ser Asn Ile Cys Cys Gly Phe Glu Ala Asn Gly Thr Val Asp Glu Ile  
 325 330 335  
 Ile Ala Ser Met Lys Leu Phe Gln Ser Asp Val Met Pro Phe Leu Lys  
 340 345 350  
 Glu Lys Gln Arg Ser Leu Leu Tyr  
 355 360

## (2) INFORMATION FOR SEQ ID NO: 7:

- (i) SEQUENCE CHARACTERISTICS:  
 (A) LENGTH: 327 amino acids  
 (B) TYPE: amino acid  
 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 7:

Met Lys Phe Gly Leu Phe Phe Leu Asn Phe Ile Asn Ser Thr Thr Val  
 1 5 10 15  
 Gln Glu Gln Ser Ile Val Arg Met Gln Glu Ile Thr Glu Tyr Val Asp  
 20 25 30  
 Lys Leu Asn Phe Glu Gln Ile Leu Val Tyr Glu Asn His Phe Ser Asp  
 35 40 45  
 Asn Gly Val Val Gly Ala Pro Leu Thr Val Ser Gly Phe Leu Leu Gly  
 50 55 60  
 Leu Thr Glu Lys Ile Lys Ile Gly Ser Leu Asn His Ile Ile Thr Thr  
 65 70 75 80  
 His His Pro Val Ala Ile Ala Glu Glu Ala Cys Leu Leu Asp Gln Leu  
 85 90 95  
 Ser Glu Gly Arg Phe Ile Leu Gly Phe Ser Asp Cys Glu Lys Lys Asp  
 100 105 110  
 Glu Met His Phe Phe Asn Arg Pro Val Glu Tyr Gln Gln Gln Leu Phe  
 115 120 125  
 Glu Glu Cys Tyr Glu Ile Ile Asn Asp Ala Leu Thr Thr Gly Tyr Cys  
 130 135 140  
 Asn Pro Asp Asn Asp Phe Tyr Ser Phe Pro Lys Ile Ser Val Asn Pro  
 145 150 155 160  
 His Ala Tyr Thr Pro Gly Gly Pro Arg Lys Tyr Val Thr Ala Thr Ser  
 165 170 175  
 His His Ile Val Glu Trp Ala Ala Lys Lys Gly Ile Pro Leu Ile Phe  
 180 185 190  
 Lys Trp Asp Asp Ser Asn Asp Val Arg Tyr Glu Tyr Ala Glu Arg Tyr  
 195 200 205

004240 29562560

Lys Ala Val Ala Asp Lys Tyr Asp Val Asp Leu Ser Glu Ile Asp His  
 210 215 220  
 Gln Leu Met Ile Leu Val Asn Tyr Asn Glu Asp Ser Asn Lys Ala Lys  
 225 230 235 240  
 Gln Glu Thr Arg Ala Phe Ile Ser Asp Tyr Val Leu Glu Met His Pro  
 245 250 255  
 Asn Glu Asn Phe Glu Asn Lys Leu Glu Glu Ile Ile Ala Glu Asn Ala  
 260 265 270  
 Val Gly Asn Tyr Thr Glu Cys Ile Thr Ala Ala Lys Leu Ala Ile Glu  
 275 280 285  
 Lys Cys Gly Ala Lys Ser Val Leu Leu Ser Phe Glu Pro Met Asn Asp  
 290 295 300  
 Leu Met Ser Gln Lys Asn Val Ile Asn Ile Val Asp Asp Asn Ile Lys  
 305 310 315 320  
 Lys Tyr His Met Glu Tyr Thr  
 325

## (2) INFORMATION FOR SEQ ID NO: 8:

- (i) SEQUENCE CHARACTERISTICS:  
 (A) LENGTH: 394 amino acids  
 (B) TYPE: amino acid  
 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

(xi) SEQUENCE DESCRIPTION: SEQ ID NO: 8:

Met Lys Gly Ile Lys Glu Tyr Asp Ser Ser Ala Ala Ile Leu Ser Asn  
 1 5 10 15  
 Ile Ile Leu Arg Ser Lys Thr Gly Met Thr Ser Tyr Val Asp Lys Gln  
 20 25 30  
 Glu Ile Thr Ala Ser Ser Glu Ile Asp Asp Leu Ile Phe Ser Ser Asp  
 35 40 45  
 Pro Leu Val Trp Ser Tyr Asp Glu Gln Glu Lys Ile Arg Lys Lys Leu  
 50 55 60  
 Val Leu Asp Ala Phe Arg Asn His Tyr Lys His Cys Arg Glu Tyr Arg  
 65 70 75 80  
 His Tyr Cys Gln Ala His Lys Val Asp Asp Asn Ile Thr Glu Ile Asp  
 85 90 95  
 Asp Ile Pro Val Phe Pro Thr Ser Val Phe Lys Phe Thr Arg Leu Leu  
 100 105 110  
 Thr Ser Gln Glu Asn Glu Ile Glu Ser Trp Phe Thr Ser Ser Gly Thr  
 115 120 125  
 Asn Gly Leu Lys Ser Gln Val Ala Arg Asp Arg Leu Ser Ile Glu Arg  
 130 135 140  
 Leu Leu Gly Ser Val Ser Tyr Gly Met Lys Tyr Val Gly Ser Trp Phe  
 145 150 155 160  
 Asp His Gln Ile Glu Leu Val Asn Leu Gly Pro Asp Arg Phe Asn Ala  
 165 170 175

004240-29662560

His Asn Ile Trp Phe Lys Tyr Val Met Ser Leu Val Glu Leu Leu Tyr  
 180 185 190  
 Pro Thr Thr Phe Thr Val Thr Glu Glu Arg Ile Asp Phe Val Lys Thr  
 195 200 205  
 Leu Asn Ser Leu Glu Arg Ile Lys Asn Gln Gly Lys Asp Leu Cys Leu  
 210 215 220  
 Ile Gly Ser Pro Tyr Phe Ile Tyr Leu Leu Cys His Tyr Met Lys Asp  
 225 230 235 240  
 Lys Lys Ile Ser Phe Ser Gly Asp Lys Ser Leu Tyr Ile Ile Thr Gly  
 245 250 255  
 Gly Gly Trp Lys Ser Tyr Glu Lys Glu Ser Leu Lys Arg Asp Asp Phe  
 260 265 270  
 Asn His Leu Leu Phe Asp Thr Phe Asn Leu Ser Asp Ile Ser Gln Ile  
 275 280 285  
 Arg Asp Ile Phe Asn Gln Val Glu Leu Asn Thr Cys Phe Phe Glu Asp  
 290 295 300  
 Glu Met Gln Arg Lys His Val Pro Pro Trp Val Tyr Ala Arg Ala Leu  
 305 310 315 320  
 Asp Pro Glu Thr Leu Lys Pro Val Pro Asp Gly Thr Pro Gly Leu Met  
 325 330 335  
 Ser Tyr Met Asp Ala Ser Ala Thr Ser Tyr Pro Ala Phe Ile Val Thr  
 340 345 350  
 Asp Asp Val Gly Ile Ile Ser Arg Glu Tyr Gly Lys Tyr Pro Gly Val  
 355 360 365  
 Leu Val Glu Ile Leu Arg Arg Val Asn Thr Arg Thr Gln Lys Gly Cys  
 370 375 380  
 Ala Leu Ser Leu Thr Glu Ala Phe Asp Ser  
 385 390

## (2) INFORMATION FOR SEQ ID NO: 9:

## (i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 3098 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: double
- (D) TOPOLOGY: circular

## (ii) MOLECULE TYPE: DNA (genomic)

## (vii) IMMEDIATE SOURCE:

- (B) CLONE: pASK75

## (viii) POSITION IN GENOME:

- (A) CHROMOSOME/SEGMENT: vector

## (ix) FEATURE:

- (A) NAME/KEY: promoter
- (B) LOCATION: 542..672
- (D) OTHER INFORMATION: /function= "beta-la promoter"  
/label= beta-la  
/citation= ([1])

00440 25866

## (ix) FEATURE:

- (A) NAME/KEY: CDS
- (B) LOCATION: 673..1530
- (D) OTHER INFORMATION: /product= "beta-la"  
/citation= ([1])

## (ix) FEATURE:

- (A) NAME/KEY: CDS
- (B) LOCATION: 1543..2163
- (D) OTHER INFORMATION: /product= "tetR"  
/citation= ([1])

## (ix) FEATURE:

- (A) NAME/KEY: misc\_feature
- (B) LOCATION: 2713..2950
- (D) OTHER INFORMATION: /function= "ORI"  
/label= ORI  
/citation= ([1])

## (ix) FEATURE:

- (A) NAME/KEY: promoter
- (B) LOCATION: 2976..3073
- (D) OTHER INFORMATION: /function= "p tetA promoter"  
/citation= ([1])

## (x) PUBLICATION INFORMATION:

- (A) AUTHORS: Skerra, A
- (B) TITLE: Use of the tetracycline promoter for the tightly regulated production of a murine antibody fragment in Escherichia coli
- (C) JOURNAL: Gene
- (D) VOLUME: 151
- (E) ISSUE: 1-2
- (F) PAGES: 131-135
- (G) DATE: 30-DEC-1994
- (K) RELEVANT RESIDUES IN SEQ ID NO: 9: FROM 1 TO 3098

## (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 9:

```
AGCTTGACCT GTGAAGTGAA AAATGGCGCA CATTGTGCGA CATTTTTTTTT GTCTGCCGTT      60
TACCGCTACT GCGTCACGGA TCTCCACGCG CCCTGTAGCG GCGCATTAAAG CGCGGCGGGT      120
GTGGTGTTA CGGCAGCGT GACCGCTACA CTTGCCAGCG CCCTAGCGCC CGCTCCTTTC      180
GCTTTCCTCC CTCCTTTCT CGCCACGTTT GCCGGCTTTC CCCGTCAAGC TCTAAATCGG      240
GGGCTCCCTT TAGGGTCCG ATTTAGTGCT TTACGGCACC TCGACCCCAA AAAACTTGAT      300
TAGGGTGATG GTTCACGTAG TGGGCCATCG CCCTGATAGA CGGTTTTTCG CCCTTTGACG      360
TTGGAGTCCA CGTCTTTTAA TAGTGGACTC TTGTTCCAAA CTGGAACAAC ACTCAACCCT      420
ATCTCGTCT ATTCTTTTGA TTTATAAGGG ATTTTGCCGA TTTCGGCCTA TTGGTTAAAA      480
AATGAGCTGA TTTAACAAAA ATTTAACGCG AATTTTAACA AAATATTAAC GCTTACAATT      540
TCAGGTGGCA CTTTTCGGGG AAATGTGCGC GGAACCCCTA TTTGTTTATT TTTCTAAATA      600
CATTCAAATA TGTATCCGCT CATGAGACAA TAACCCTGAT AAATGCTTCA ATAATATTGA      660
AAAAGGAAGA GT ATG AGT ATT CAA CAT TTC CGT GTC GCC CTT ATT CCC      708
          Met Ser Ile Gln His Phe Arg Val Ala Leu Ile Pro
          395                400                405
```

004240 29662560

TTT TTT GCG GCA TTT TGC CTT CCT GTT TTT GCT CAC CCA GAA ACG CTG 756  
 Phe Phe Ala Ala Phe Cys Leu Pro Val Phe Ala His Pro Glu Thr Leu  
 410 415 420  
 GTG AAA GTA AAA GAT GCT GAA GAT CAG TTG GGT GCA CGA GTG GGT TAC 804  
 Val Lys Val Lys Asp Ala Glu Asp Gln Leu Gly Ala Arg Val Gly Tyr  
 425 430 435  
 ATC GAA CTG GAT CTC AAC AGC GGT AAG ATC CTT GAG AGT TTT CGC CCC 852  
 Ile Glu Leu Asp Leu Asn Ser Gly Lys Ile Leu Glu Ser Phe Arg Pro  
 440 445 450  
 GAA GAA CGT TTT CCA ATG ATG AGC ACT TTT AAA GTT CTG CTA TGT GGC 900  
 Glu Glu Arg Phe Pro Met Met Ser Thr Phe Lys Val Leu Leu Cys Gly  
 455 460 465 470  
 GCG GTA TTA TCC CGT ATT GAC GCC GGG CAA GAG CAA CTC GGT CGC CGC 948  
 Ala Val Leu Ser Arg Ile Asp Ala Gly Gln Glu Gln Leu Gly Arg Arg  
 475 480 485  
 ATA CAC TAT TCT CAG AAT GAC TTG GTT GAG TAC TCA CCA GTC ACA GAA 996  
 Ile His Tyr Ser Gln Asn Asp Leu Val Glu Tyr Ser Pro Val Thr Glu  
 490 495 500  
 AAG CAT CTT ACG GAT GGC ATG ACA GTA AGA GAA TTA TGC AGT GCT GCC 1044  
 Lys His Leu Thr Asp Gly Met Thr Val Arg Glu Leu Cys Ser Ala Ala  
 505 510 515  
 ATA ACC ATG AGT GAT AAC ACT GCG GCC AAC TTA CTT CTG ACA ACG ATC 1092  
 Ile Thr Met Ser Asp Asn Thr Ala Ala Asn Leu Leu Thr Thr Ile  
 520 525 530  
 GGA GGA CCG AAG GAG CTA ACC GCT TTT TTG CAC AAC ATG GGG GAT CAT 1140  
 Gly Gly Pro Lys Glu Leu Thr Ala Phe Leu His Asn Met Gly Asp His  
 535 540 545 550  
 GTA ACT CGC CTT GAT CGT TGG GAA CCG GAG CTG AAT GAA GCC ATA CCA 1188  
 Val Thr Arg Leu Asp Arg Trp Glu Pro Glu Leu Asn Glu Ala Ile Pro  
 555 560 565  
 AAC GAC GAG CGT GAC ACC ACG ATG CCT GTA GCA ATG GCA ACA ACG TTG 1236  
 Asn Asp Glu Arg Asp Thr Thr Met Pro Val Ala Met Ala Thr Thr Leu  
 570 575 580  
 CGC AAA CTA TTA ACT GGC GAA CTA CTT ACT CTA GCT TCC CGG CAA CAA 1284  
 Arg Lys Leu Leu Thr Gly Glu Leu Leu Thr Leu Ala Ser Arg Gln Gln  
 585 590 595  
 TTG ATA GAC TGG ATG GAG GCG GAT AAA GTT GCA GGA CCA CTT CTG CGC 1332  
 Leu Ile Asp Trp Met Glu Ala Asp Lys Val Ala Gly Pro Leu Leu Arg  
 600 605 610  
 TCG GCC CTT CCG GCT GGC TGG TTT ATT GCT GAT AAA TCT GGA GCC GGT 1380  
 Ser Ala Leu Pro Ala Gly Trp Phe Ile Ala Asp Lys Ser Gly Ala Gly  
 615 620 625 630  
 GAG CGT GGC TCT CGC GGT ATC ATT GCA GCA CTG GGG CCA GAT GGT AAG 1428  
 Glu Arg Gly Ser Arg Gly Ile Ile Ala Ala Leu Gly Pro Asp Gly Lys  
 635 640 645  
 CCC TCC CGT ATC GTA GTT ATC TAC ACG ACG GGG AGT CAG GCA ACT ATG 1476  
 Pro Ser Arg Ile Val Val Ile Tyr Thr Thr Gly Ser Gln Ala Thr Met  
 650 655 660  
 GAT GAA CGA AAT AGA CAG ATC GCT GAG ATA GGT GCC TCA CTG ATT AAG 1524  
 Asp Glu Arg Asn Arg Gln Ile Ala Glu Ile Gly Ala Ser Leu Ile Lys  
 665 670 675

004240" 29662560

CAT TGG TAGGAATTAA TG ATG TCT CGT TTA GAT AAA AGT AAA GTG ATT 1572  
 His Trp Met Ser Arg Leu Asp Lys Ser Lys Val Ile  
 680 1 5 10

AAC AGC GCA TTA GAG CTG CTT AAT GAG GTC GGA ATC GAA GGT TTA ACA 1620  
 Asn Ser Ala Leu Glu Leu Leu Asn Glu Val Gly Ile Glu Gly Leu Thr  
 15 20 25

ACC CGT AAA CTC GCC CAG AAG CTA GGT GTA GAG CAG CCT ACA TTG TAT 1668  
 Thr Arg Lys Leu Ala Gln Lys Leu Gly Val Glu Gln Pro Thr Leu Tyr  
 30 35 40

TGG CAT GTA AAA AAT AAG CGG GCT TTG CTC GAC GCC TTA GCC ATT GAG 1716  
 Trp His Val Lys Asn Lys Arg Ala Leu Leu Asp Ala Leu Ala Ile Glu  
 45 50 55

ATG TTA GAT AGG CAC CAT ACT CAC TTT TGC CCT TTA GAA GGG GAA AGC 1764  
 Met Leu Asp Arg His His Thr His Phe Cys Pro Leu Glu Gly Glu Ser  
 60 65 70

TGG CAA GAT TTT TTA CGT AAT AAC GCT AAA AGT TTT AGA TGT GCT TTA 1812  
 Trp Gln Asp Phe Leu Arg Asn Asn Ala Lys Ser Phe Arg Cys Ala Leu  
 75 80 85 90

CTA AGT CAT CGC GAT GGA GCA AAA GTA CAT TTA GGT ACA CGG CCT ACA 1860  
 Leu Ser His Arg Asp Gly Ala Lys Val His Leu Gly Thr Arg Pro Thr  
 95 100 105

GAA AAA CAG TAT GAA ACT CTC GAA AAT CAA TTA GCC TTT TTA TGC CAA 1908  
 Glu Lys Gln Tyr Glu Thr Leu Glu Asn Gln Leu Ala Phe Leu Cys Gln  
 110 115 120

CAA GGT TTT TCA CTA GAG AAT GCA TTA TAT GCA CTC AGC GCA GTG GGG 1956  
 Gln Gly Phe Ser Leu Glu Asn Ala Leu Tyr Ala Leu Ser Ala Val Gly  
 125 130 135

CAT TTT ACT TTA GGT TGC GTA TTG GAA GAT CAA GAG CAT CAA GTC GCT 2004  
 His Phe Thr Leu Gly Cys Val Leu Glu Asp Gln Glu His Gln Val Ala  
 140 145 150

AAA GAA GAA AGG GAA ACA CCT ACT ACT GAT AGT ATG CCG CCA TTA TTA 2052  
 Lys Glu Glu Arg Glu Thr Pro Thr Thr Asp Ser Met Pro Pro Leu Leu  
 155 160 165 170

CGA CAA GCT ATC GAA TTA TTT GAT CAC CAA GGT GCA GAG CCA GCC TTC 2100  
 Arg Gln Ala Ile Glu Leu Phe Asp His Gln Gly Ala Glu Pro Ala Phe  
 175 180 185

TTA TTC GGC CTT GAA TTG ATC ATA TGC GGA TTA GAA AAA CAA CTT AAA 2148  
 Leu Phe Gly Leu Glu Leu Ile Ile Cys Gly Leu Glu Lys Gln Leu Lys  
 190 195 200

TGT GAA AGT GGG TCT TAAAGCAGC ATAACCTTTT TCCGTGATGG TAACTTCACT 2203  
 Cys Glu Ser Gly Ser  
 205

AGTTTAAAAG GATCTAGGTG AAGATCCTTT TTGATAATCT CATGACCAAA ATCCCTTAAC 2263

GTGAGTTTTT GTTCCACTGA GCGTCAGACC CCGTAGAAAA GATCAAAGGA TCTTCTTGAG 2323

ATCCTTTTTT TCTGCGCGTA ATCTGCTGCT TGCAAACAAA AAAACCACCG CTACCAGCGG 2383

TGGTTTGTTT GCCGGATCAA GAGCTACCAA CTCTTTTTCC GAAGGTAACCT GGCTTCAGCA 2443

GAGCGCAGAT ACCAAATACT GTCCTTCTAG TGTAGCCGTA GTTAGGCCAC CACTTCAAGA 2503

ACTCTGTAGC ACCGCCTACA TACCTCGCTC TGCTAATCCT GTTACCAGTG GCTGCTGCCA 2563

004240 29662560

GTGGCGATAA GTCGTGTCTT ACCGGGTTGG ACTCAAGACG ATAGTTACCG GATAAGGCGC 2623  
 AGCGGTCGGG CTGAACGGGG GGTTCGTGCA CACAGCCCAG CTTGGAGCGA ACGACCTACA 2683  
 CCGAACTGAG ATACCTACAG CGTGAGCTAT GAGAAAGCGC CACGCTTCCC GAAGGGAGAA 2743  
 AGGCGGACAG GTATCCGGTA AGCGGCAGGG TCGGAACAGG AGAGCGCACG AGGGAGCTTC 2803  
 CAGGGGGAAA CGCCTGGTAT CTTTATAGTC CTGTCCGGTT TCGCCACCTC TGACTTGAGC 2863  
 GTCGATTTTT GTGATGCTCG TCAGGGGGGC GGAGCCTATG GAAAAAGGCC AGCAACGCGG 2923  
 CCTTTTTACG GTTCCTGGCC TTTTGCTGGC CTTTGTCTCA CATGACCCGA CACCATCGAA 2983  
 TGGCCAGATG ATTAATTCCT AATTTTGTGTT GACACTCTAT CATTGATAGA GTTATTTTAC 3043  
 CACTCCCTAT CAGTGATAGA GAAAAGTGAA ATGAATAGTT CGACAAAAAT CTAGA 3098

## (2) INFORMATION FOR SEQ ID NO: 10:

- (i) SEQUENCE CHARACTERISTICS:  
 (A) LENGTH: 286 amino acids  
 (B) TYPE: amino acid  
 (D) TOPOLOGY: linear

- (ii) MOLECULE TYPE: protein  
 (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 10:

Met Ser Ile Gln His Phe Arg Val Ala Leu Ile Pro Phe Phe Ala Ala  
 1 5 10 15  
 Phe Cys Leu Pro Val Phe Ala His Pro Glu Thr Leu Val Lys Val Lys  
 20 25 30  
 Asp Ala Glu Asp Gln Leu Gly Ala Arg Val Gly Tyr Ile Glu Leu Asp  
 35 40 45  
 Leu Asn Ser Gly Lys Ile Leu Glu Ser Phe Arg Pro Glu Glu Arg Phe  
 50 55 60  
 Pro Met Met Ser Thr Phe Lys Val Leu Leu Cys Gly Ala Val Leu Ser  
 65 70 75 80  
 Arg Ile Asp Ala Gly Gln Glu Gln Leu Gly Arg Arg Ile His Tyr Ser  
 85 90 95  
 Gln Asn Asp Leu Val Glu Tyr Ser Pro Val Thr Glu Lys His Leu Thr  
 100 105 110  
 Asp Gly Met Thr Val Arg Glu Leu Cys Ser Ala Ala Ile Thr Met Ser  
 115 120 125  
 Asp Asn Thr Ala Ala Asn Leu Leu Leu Thr Thr Ile Gly Gly Pro Lys  
 130 135 140  
 Glu Leu Thr Ala Phe Leu His Asn Met Gly Asp His Val Thr Arg Leu  
 145 150 155 160  
 Asp Arg Trp Glu Pro Glu Leu Asn Glu Ala Ile Pro Asn Asp Glu Arg  
 165 170 175  
 Asp Thr Thr Met Pro Val Ala Met Ala Thr Thr Leu Arg Lys Leu Leu  
 180 185 190  
 Thr Gly Glu Leu Leu Thr Leu Ala Ser Arg Gln Gln Leu Ile Asp Trp  
 195 200 205

004240 29662560

Met Glu Ala Asp Lys Val Ala Gly Pro Leu Leu Arg Ser Ala Leu Pro  
 210 215 220  
 Ala Gly Trp Phe Ile Ala Asp Lys Ser Gly Ala Gly Glu Arg Gly Ser  
 225 230 235 240  
 Arg Gly Ile Ile Ala Ala Leu Gly Pro Asp Gly Lys Pro Ser Arg Ile  
 245 250 255  
 Val Val Ile Tyr Thr Thr Gly Ser Gln Ala Thr Met Asp Glu Arg Asn  
 260 265 270  
 Arg Gln Ile Ala Glu Ile Gly Ala Ser Leu Ile Lys His Trp  
 275 280 285

## (2) INFORMATION FOR SEQ ID NO: 11:

## (i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 207 amino acids  
 (B) TYPE: amino acid  
 (D) TOPOLOGY: linear

## (ii) MOLECULE TYPE: protein

## (xi) SEQUENCE DESCRIPTION: SEQ ID NO: 11:

Met Ser Arg Leu Asp Lys Ser Lys Val Ile Asn Ser Ala Leu Glu Leu  
 1 5 10 15  
 Leu Asn Glu Val Gly Ile Glu Gly Leu Thr Thr Arg Lys Leu Ala Gln  
 20 25 30  
 Lys Leu Gly Val Glu Gln Pro Thr Leu Tyr Trp His Val Lys Asn Lys  
 35 40 45  
 Arg Ala Leu Leu Asp Ala Leu Ala Ile Glu Met Leu Asp Arg His His  
 50 55 60  
 Thr His Phe Cys Pro Leu Glu Gly Glu Ser Trp Gln Asp Phe Leu Arg  
 65 70 75 80  
 Asn Asn Ala Lys Ser Phe Arg Cys Ala Leu Leu Ser His Arg Asp Gly  
 85 90 95  
 Ala Lys Val His Leu Gly Thr Arg Pro Thr Glu Lys Gln Tyr Glu Thr  
 100 105 110  
 Leu Glu Asn Gln Leu Ala Phe Leu Cys Gln Gln Gly Phe Ser Leu Glu  
 115 120 125  
 Asn Ala Leu Tyr Ala Leu Ser Ala Val Gly His Phe Thr Leu Gly Cys  
 130 135 140  
 Val Leu Glu Asp Gln Glu His Gln Val Ala Lys Glu Glu Arg Glu Thr  
 145 150 155 160  
 Pro Thr Thr Asp Ser Met Pro Pro Leu Leu Arg Gln Ala Ile Glu Leu  
 165 170 175  
 Phe Asp His Gln Gly Ala Glu Pro Ala Phe Leu Phe Gly Leu Glu Leu  
 180 185 190  
 Ile Ile Cys Gly Leu Glu Lys Gln Leu Lys Cys Glu Ser Gly Ser  
 195 200 205

004240 29523560